

PC12 cells with a construct containing an EGR promoter linked to SEAP reporter, activation of PC12 cells by polypeptide of the present invention can be assessed.

[01168] The EGR/SEAP reporter construct can be assembled by the following protocol. The EGR-1 promoter sequence (-633 to +1)(Sakamoto K et al., Oncogene 6:867-871 (1991)) can be PCR amplified from human genomic DNA using the following primers:

5' GCGCTCGAGGGATGACAGCGATAGAACCCCGG -3' (SEQ ID NO: 6)

5' GCGAAGCTTCGCGACTCCCCGGATCCGCCTC-3' (SEQ ID NO: 7)

[01169] Using the GAS:SEAP/Neo vector produced in Example 31, EGR1 amplified product can then be inserted into this vector. Linearize the GAS:SEAP/Neo vector using restriction enzymes XhoI/HindIII, removing the GAS/SV40 stuffer. Restrict the EGR1 amplified product with these same enzymes. Ligate the vector and the EGR1 promoter.

[01170] To prepare 96 well-plates for cell culture, two mls of a coating solution (1:30 dilution of collagen type I (Upstate Biotech Inc. Cat#08-115) in 30% ethanol (filter sterilized)) is added per one 10 cm plate or 50 ml per well of the 96-well plate, and allowed to air dry for 2 hr.

[01171] PC12 cells are routinely grown in RPMI-1640 medium (Bio Whittaker) containing 10% horse serum (JRH BIOSCIENCES, Cat. # 12449-78P), 5% heat-inactivated fetal bovine serum (FBS) supplemented with 100 units/ml penicillin and 100 ug/ml streptomycin on a precoated 10 cm tissue culture dish. One to four split is done every three to four days. Cells are removed from the plates by scraping and resuspended with pipetting up and down for more than 15 times.

[01172] Transfect the EGR/SEAP/Neo construct into PC12 using the Lipofectamine protocol described in Example 30. EGR-SEAP/PC12 stable cells are obtained by growing the cells in 300 ug/ml G418. The G418-free medium is used for routine growth but every one to two months, the cells should be re-grown in 300 ug/ml G418 for couple of passages.

[01173] To assay for neuronal activity, a 10 cm plate with cells around 70 to 80% confluent is screened by removing the old medium. Wash the cells once with PBS (Phosphate buffered saline). Then starve the cells in low serum medium (RPMI-1640 containing 1% horse serum and 0.5% FBS with antibiotics) overnight.

[01174] The next morning, remove the medium and wash the cells with PBS. Scrape off the cells from the plate, suspend the cells well in 2 ml low serum medium. Count the cell number and add more low serum medium to reach final cell density as  $5 \times 10^5$  cells/ml.

[01175] Add 200  $\mu$ l of the cell suspension to each well of 96-well plate (equivalent to  $1 \times 10^5$  cells/well). Add 50  $\mu$ l supernatant produced by Example 30, 37 degree C for 48 to 72 hr. As a positive control, a growth factor known to activate PC12 cells through EGR can be used, such as 50 ng/ $\mu$ l of Neuronal Growth Factor (NGF). Over fifty-fold induction of SEAP is typically seen in the positive control wells. SEAP assay the supernatant according to Example 36.

### ***Example 35: High-Throughput Screening Assay for T-cell Activity***

[01176] NF-KB (Nuclear Factor KB) is a transcription factor activated by a wide variety of agents including the inflammatory cytokines IL-1 and TNF, CD30 and CD40, lymphotoxin-alpha and lymphotoxin-beta, by exposure to LPS or thrombin, and by expression of certain viral gene products. As a transcription factor, NF-KB regulates the expression of genes involved in immune cell activation, control of apoptosis (NF- KB appears to shield cells from apoptosis), B and T-cell development, anti-viral and antimicrobial responses, and multiple stress responses.

[01177] In non-stimulated conditions, NF- KB is retained in the cytoplasm with I-KB (Inhibitor KB). However, upon stimulation, I- KB is phosphorylated and degraded, causing NF- KB to shuttle to the nucleus, thereby activating transcription of target genes. Target genes activated by NF- KB include IL-2, IL-6, GM-CSF, ICAM-1 and class 1 MHC.

[01178] Due to its central role and ability to respond to a range of stimuli, reporter constructs utilizing the NF-KB promoter element are used to screen the supernatants produced in Example 30. Activators or inhibitors of NF-KB would be useful in treating, preventing, and/or diagnosing diseases. For example, inhibitors of NF-KB could be used to treat those diseases related to the acute or chronic activation of NF-KB, such as rheumatoid arthritis.

[01179] To construct a vector containing the NF-KB promoter element, a PCR based strategy is employed. The upstream primer contains four tandem copies of the NF-KB binding site (GGGGACTTTCCC) (SEQ ID NO: 8), 18 bp of sequence complementary to the 5' end of the SV40 early promoter sequence, and is flanked with an XhoI site:

5':GCGGCCTCGAGGGGACTTTCCCGGGGACTTTCCGGGGACTTTCCGGGGACT  
TTCCATCCTGCCATCTCAATTAG:3' (SEQ ID NO: 9)

[01180] The downstream primer is complementary to the 3' end of the SV40 promoter and is flanked with a Hind III site:

5':GCGGCAAGCTTTTTGCAAAGCCTAGGC:3' (SEQ ID NO: 4)

[01181] PCR amplification is performed using the SV40 promoter template present in the pB-gal:promoter plasmid obtained from Clontech. The resulting PCR fragment is digested with XhoI and Hind III and subcloned into BLSK2-. (Stratagene) Sequencing with the T7 and T3 primers confirms the insert contains the following sequence:

5':CTCGAGGGGACTTTCCCGGGGACTTTCCGGGGACTTTCCGGGGACTTTCCAT  
CTGCCATCTCAATTAGTCAGCAACCATAGTCCCGCCCCCTAACTCCGCCCCATCC  
CGCCCCCTAACTCCGCCCCAGTTCCGCCCCATTCTCCGCCCCATGGCTGACTAATT  
TTTTTTATTTATGCAGAGGCCGAGGCCGCTCGGCCTCTGAGCTATTCCAGAA  
GTAGTGAGGAGGCTTTTTTGGAGGCCTAGGCTTTTGCAAAAAGCTT:3' (SEQ  
ID NO: 10)

[01182] Next, replace the SV40 minimal promoter element present in the pSEAP2-promoter plasmid (Clontech) with this NF-KB/SV40 fragment using XhoI and HindIII. However, this vector does not contain a neomycin resistance gene, and therefore, is not preferred for mammalian expression systems.

[01183] In order to generate stable mammalian cell lines, the NF-KB/SV40/SEAP cassette is removed from the above NF-KB/SEAP vector using restriction enzymes SalI and NotI, and inserted into a vector containing neomycin resistance. Particularly, the NF-KB/SV40/SEAP cassette was inserted into pGFP-1 (Clontech), replacing the GFP gene, after restricting pGFP-1 with SalI and NotI.

[01184] Once NF-KB/SV40/SEAP/Neo vector is created, stable Jurkat T-cells are created and maintained according to the protocol described in Example 32. Similarly, the method for assaying supernatants with these stable Jurkat T-cells is also described in

Example 32. As a positive control, exogenous TNF alpha (0.1, 1, 10 ng) is added to wells H9, H10, and H11, with a 5-10 fold activation typically observed.

### ***Example 36: Assay for SEAP Activity***

[01185] As a reporter molecule for the assays described in Examples 32-35, SEAP activity is assayed using the Tropix Phospho-light Kit (Cat. BP-400) according to the following general procedure. The Tropix Phospho-light Kit supplies the Dilution, Assay, and Reaction Buffers used below.

[01186] Prime a dispenser with the 2.5x Dilution Buffer and dispense 15 ul of 2.5x dilution buffer into Optiplates containing 35 ul of a supernatant. Seal the plates with a plastic sealer and incubate at 65 degree C for 30 min. Separate the Optiplates to avoid uneven heating.

[01187] Cool the samples to room temperature for 15 minutes. Empty the dispenser and prime with the Assay Buffer. Add 50 ml Assay Buffer and incubate at room temperature 5 min. Empty the dispenser and prime with the Reaction Buffer (see the Table below). Add 50 ul Reaction Buffer and incubate at room temperature for 20 minutes. Since the intensity of the chemiluminescent signal is time dependent, and it takes about 10 minutes to read 5 plates on a luminometer, thus one should treat 5 plates at each time and start the second set 10 minutes later.

[01188] Read the relative light unit in the luminometer. Set H12 as blank, and print the results. An increase in chemiluminescence indicates reporter activity.

#### **Reaction Buffer Formulation:**

# of plates	Rxn buffer diluent (ml)	CSPD (ml)
10	60	3
11	65	3.25
12	70	3.5
13	75	3.75
14	80	4
15	85	4.25
16	90	4.5
17	95	4.75

18	100	5
19	105	5.25
20	110	5.5
21	115	5.75
22	120	6
23	125	6.25
24	130	6.5
25	135	6.75
26	140	7
27	145	7.25
28	150	7.5
29	155	7.75
30	160	8
31	165	8.25
32	170	8.5
33	175	8.75
34	180	9
35	185	9.25
36	190	9.5
37	195	9.75
38	200	10
39	205	10.25
40	210	10.5
41	215	10.75
42	220	11
43	225	11.25
44	230	11.5
45	235	11.75
46	240	12
47	245	12.25
48	250	12.5
49	255	12.75
50	260	13

---

***Example 37: High-Throughput Screening Assay Identifying Changes in  
Small Molecule Concentration and Membrane Permeability***

[01189] Binding of a ligand to a receptor is known to alter intracellular levels of small molecules, such as calcium, potassium, sodium, and pH, as well as alter membrane potential. These alterations can be measured in an assay to identify supernatants which bind to receptors of a particular cell. Although the following protocol describes an assay for calcium, this protocol can easily be modified to detect changes in potassium, sodium, pH, membrane potential, or any other small molecule which is detectable by a fluorescent probe.

[01190] The following assay uses Fluorometric Imaging Plate Reader ("FLIPR") to measure changes in fluorescent molecules (Molecular Probes) that bind small molecules. Clearly, any fluorescent molecule detecting a small molecule can be used instead of the calcium fluorescent molecule, fluo-4 (Molecular Probes, Inc.; catalog no. F-14202), used here.

[01191] For adherent cells, seed the cells at 10,000 -20,000 cells/well in a Co-star black 96-well plate with clear bottom. The plate is incubated in a CO<sub>2</sub> incubator for 20 hours. The adherent cells are washed two times in Biotek washer with 200 ul of HBSS (Hank's Balanced Salt Solution) leaving 100 ul of buffer after the final wash.

[01192] A stock solution of 1 mg/ml fluo-4 is made in 10% pluronic acid DMSO. To load the cells with fluo-4, 50 ul of 12 ug/ml fluo-4 is added to each well. The plate is incubated at 37 degrees C in a CO<sub>2</sub> incubator for 60 min. The plate is washed four times in the Biotek washer with HBSS leaving 100 ul of buffer.

[01193] For non-adherent cells, the cells are spun down from culture media. Cells are re-suspended to  $2-5 \times 10^6$  cells/ml with HBSS in a 50-ml conical tube. 4 ul of 1 mg/ml fluo-4 solution in 10% pluronic acid DMSO is added to each ml of cell suspension. The tube is then placed in a 37 degrees C water bath for 30-60 min. The cells are washed twice with HBSS, resuspended to  $1 \times 10^6$  cells/ml, and dispensed into a microplate, 100 ul/well. The plate is centrifuged at 1000 rpm for 5 min. The plate is then washed once in Denley Cell Wash with 200 ul, followed by an aspiration step to 100 ul final volume.

[01194] For a non-cell based assay, each well contains a fluorescent molecule, such as fluo-4. The supernatant is added to the well, and a change in fluorescence is detected.

[01195] To measure the fluorescence of intracellular calcium, the FLIPR is set for the following parameters: (1) System gain is 300-800 mW; (2) Exposure time is 0.4 second; (3) Camera F/stop is F/2; (4) Excitation is 488 nm; (5) Emission is 530 nm; and

(6) Sample addition is 50 ul. Increased emission at 530 nm indicates an extracellular signaling event caused by the a molecule, either polypeptide of the present invention or a molecule induced by polypeptide of the present invention, which has resulted in an increase in the intracellular  $\text{Ca}^{++}$  concentration.

***Example 38: High-Throughput Screening Assay Identifying Tyrosine Kinase Activity***

[01196] The Protein Tyrosine Kinases (PTK) represent a diverse group of transmembrane and cytoplasmic kinases. Within the Receptor Protein Tyrosine Kinase (RPTK) group are receptors for a range of mitogenic and metabolic growth factors including the PDGF, FGF, EGF, NGF, HGF and Insulin receptor subfamilies. In addition there are a large family of RPTKs for which the corresponding ligand is unknown. Ligands for RPTKs include mainly secreted small proteins, but also membrane-bound and extracellular matrix proteins.

[01197] Activation of RPTK by ligands involves ligand-mediated receptor dimerization, resulting in transphosphorylation of the receptor subunits and activation of the cytoplasmic tyrosine kinases. The cytoplasmic tyrosine kinases include receptor associated tyrosine kinases of the src-family (e.g., src, yes, lck, lyn, fyn) and non-receptor linked and cytosolic protein tyrosine kinases, such as the Jak family, members of which mediate signal transduction triggered by the cytokine superfamily of receptors (e.g., the Interleukins, Interferons, GM-CSF, and Leptin).

[01198] Because of the wide range of known factors capable of stimulating tyrosine kinase activity, identifying whether polypeptide of the present invention or a molecule induced by polypeptide of the present invention is capable of activating tyrosine kinase signal transduction pathways is of interest. Therefore, the following protocol is designed to identify such molecules capable of activating the tyrosine kinase signal transduction pathways.

[01199] Seed target cells (e.g., primary keratinocytes) at a density of approximately 25,000 cells per well in a 96 well Loprodyne Silent Screen Plates purchased from Nalge Nunc (Naperville, IL). The plates are sterilized with two 30 minute rinses with 100% ethanol, rinsed with water and dried overnight. Some plates are coated for 2 hr with 100

ml of cell culture grade type I collagen (50 mg/ml), gelatin (2%) or polylysine (50 mg/ml), all of which can be purchased from Sigma Chemicals (St. Louis, MO) or 10% Matrigel purchased from Becton Dickinson (Bedford, MA), or calf serum, rinsed with PBS and stored at 4 degree C. Cell growth on these plates is assayed by seeding 5,000 cells/well in growth medium and indirect quantitation of cell number through use of alamarBlue as described by the manufacturer Alamar Biosciences, Inc. (Sacramento, CA) after 48 hr. Falcon plate covers #3071 from Becton Dickinson (Bedford, MA) are used to cover the Loprodyne Silent Screen Plates. Falcon Microtest III cell culture plates can also be used in some proliferation experiments.

**[01200]** To prepare extracts, A431 cells are seeded onto the nylon membranes of Loprodyne plates (20,000/200ml/well) and cultured overnight in complete medium. Cells are quiesced by incubation in serum-free basal medium for 24 hr. After 5-20 minutes treatment with EGF (60ng/ml) or 50 ul of the supernatant produced in Example 30, the medium was removed and 100 ml of extraction buffer ((20 mM HEPES pH 7.5, 0.15 M NaCl, 1% Triton X-100, 0.1% SDS, 2 mM Na<sub>3</sub>VO<sub>4</sub>, 2 mM Na<sub>4</sub>P<sub>2</sub>O<sub>7</sub> and a cocktail of protease inhibitors (# 1836170) obtained from Boehringer Mannheim (Indianapolis, IN)) is added to each well and the plate is shaken on a rotating shaker for 5 minutes at 4°C. The plate is then placed in a vacuum transfer manifold and the extract filtered through the 0.45 mm membrane bottoms of each well using house vacuum. Extracts are collected in a 96-well catch/assay plate in the bottom of the vacuum manifold and immediately placed on ice. To obtain extracts clarified by centrifugation, the content of each well, after detergent solubilization for 5 minutes, is removed and centrifuged for 15 minutes at 4 degree C at 16,000 x g.

**[01201]** Test the filtered extracts for levels of tyrosine kinase activity. Although many methods of detecting tyrosine kinase activity are known, one method is described here.

**[01202]** Generally, the tyrosine kinase activity of a supernatant is evaluated by determining its ability to phosphorylate a tyrosine residue on a specific substrate (a biotinylated peptide). Biotinylated peptides that can be used for this purpose include PSK1 (corresponding to amino acids 6-20 of the cell division kinase cdc2-p34) and PSK2 (corresponding to amino acids 1-17 of gastrin). Both peptides are substrates for a range of tyrosine kinases and are available from Boehringer Mannheim.



[01203] The tyrosine kinase reaction is set up by adding the following components in order. First, add 10ul of 5uM Biotinylated Peptide, then 10ul ATP/Mg<sub>2</sub><sup>+</sup> (5mM ATP/50mM MgCl<sub>2</sub>), then 10ul of 5x Assay Buffer (40mM imidazole hydrochloride, pH7.3, 40 mM beta-glycerophosphate, 1mM EGTA, 100mM MgCl<sub>2</sub>, 5 mM MnCl<sub>2</sub>, 0.5 mg/ml BSA), then 5ul of Sodium Vanadate (1mM), and then 5ul of water. Mix the components gently and preincubate the reaction mix at 30 degree C for 2 min. Initiate the reaction by adding 10ul of the control enzyme or the filtered supernatant.

[01204] The tyrosine kinase assay reaction is then terminated by adding 10 ul of 120mM EDTA and place the reactions on ice.

[01205] Tyrosine kinase activity is determined by transferring 50 ul aliquot of reaction mixture to a microtiter plate (MTP) module and incubating at 37 degree C for 20 min. This allows the streptavidin coated 96 well plate to associate with the biotinylated peptide. Wash the MTP module with 300ul/well of PBS four times. Next add 75 ul of anti-phosphotyrosine antibody conjugated to horse radish peroxidase (anti-P-Tyr-POD(0.5u/ml)) to each well and incubate at 37 degree C for one hour. Wash the well as above.

[01206] Next add 100ul of peroxidase substrate solution (Boehringer Mannheim) and incubate at room temperature for at least 5 mins (up to 30 min). Measure the absorbance of the sample at 405 nm by using ELISA reader. The level of bound peroxidase activity is quantitated using an ELISA reader and reflects the level of tyrosine kinase activity.

### ***Example 39: High-Throughput Screening Assay Identifying Phosphorylation Activity***

[01207] As a potential alternative and/or complement to the assay of protein tyrosine kinase activity described in Example 38, an assay which detects activation (phosphorylation) of major intracellular signal transduction intermediates can also be used. For example, as described below one particular assay can detect tyrosine phosphorylation of the Erk-1 and Erk-2 kinases. However, phosphorylation of other molecules, such as Raf, JNK, p38 MAP, Map kinase kinase (MEK), MEK kinase, Src, Muscle specific kinase (MuSK), IRAK, Tec, and Janus, as well as any other

phosphoserine, phosphotyrosine, or phosphothreonine molecule, can be detected by substituting these molecules for Erk-1 or Erk-2 in the following assay.

[01208] Specifically, assay plates are made by coating the wells of a 96-well ELISA plate with 0.1ml of protein G (1ug/ml) for 2 hr at room temp, (RT). The plates are then rinsed with PBS and blocked with 3% BSA/PBS for 1 hr at RT. The protein G plates are then treated with 2 commercial monoclonal antibodies (100ng/well) against Erk-1 and Erk-2 (1 hr at RT) (Santa Cruz Biotechnology). (To detect other molecules, this step can easily be modified by substituting a monoclonal antibody detecting any of the above described molecules.) After 3-5 rinses with PBS, the plates are stored at 4 degree C until use.

[01209] A431 cells are seeded at 20,000/well in a 96-well Loprodyne filterplate and cultured overnight in growth medium. The cells are then starved for 48 hr in basal medium (DMEM) and then treated with EGF (6ng/well) or 50 ul of the supernatants obtained in Example 30 for 5-20 minutes. The cells are then solubilized and extracts filtered directly into the assay plate.

[01210] After incubation with the extract for 1 hr at RT, the wells are again rinsed. As a positive control, a commercial preparation of MAP kinase (10ng/well) is used in place of A431 extract. Plates are then treated with a commercial polyclonal (rabbit) antibody (1ug/ml) which specifically recognizes the phosphorylated epitope of the Erk-1 and Erk-2 kinases (1 hr at RT). This antibody is biotinylated by standard procedures. The bound polyclonal antibody is then quantitated by successive incubations with Europium-streptavidin and Europium fluorescence enhancing reagent in the Wallac DELFIA instrument (time-resolved fluorescence). An increased fluorescent signal over background indicates a phosphorylation by polypeptide of the present invention or a molecule induced by polypeptide of the present invention.

#### ***Example 40: Assay for the Stimulation of Bone Marrow CD34+ Cell Proliferation***

[01211] This assay is based on the ability of human CD34+ to proliferate in the presence of hematopoietic growth factors and evaluates the ability of isolated polypeptides expressed in mammalian cells to stimulate proliferation of CD34+ cells.

[01212] It has been previously shown that most mature precursors will respond to only a single signal. More immature precursors require at least two signals to respond. Therefore, to test the effect of polypeptides on hematopoietic activity of a wide range of progenitor cells, the assay contains a given polypeptide in the presence or absence of other hematopoietic growth factors. Isolated cells are cultured for 5 days in the presence of Stem Cell Factor (SCF) in combination with tested sample. SCF alone has a very limited effect on the proliferation of bone marrow (BM) cells, acting in such conditions only as a "survival" factor. However, combined with any factor exhibiting stimulatory effect on these cells (e.g., IL-3), SCF will cause a synergistic effect. Therefore, if the tested polypeptide has a stimulatory effect on hematopoietic progenitors, such activity can be easily detected. Since normal BM cells have a low level of cycling cells, it is likely that any inhibitory effect of a given polypeptide, or agonists or antagonists thereof, might not be detected. Accordingly, assays for an inhibitory effect on progenitors is preferably tested in cells that are first subjected to *in vitro* stimulation with SCF+IL-3, and then contacted with the compound that is being evaluated for inhibition of such induced proliferation.

[01213] Briefly, CD34+ cells are isolated using methods known in the art. The cells are thawed and resuspended in medium (QBSF 60 serum-free medium with 1% L-glutamine (500ml) Quality Biological, Inc., Gaithersburg, MD Cat# 160-204-101). After several gentle centrifugation steps at 200 x g, cells are allowed to rest for one hour. The cell count is adjusted to  $2.5 \times 10^5$  cells/ml. During this time, 100  $\mu$ l of sterile water is added to the peripheral wells of a 96-well plate. The cytokines that can be tested with a given polypeptide in this assay is rhSCF (R&D Systems, Minneapolis, MN, Cat# 255-SC) at 50 ng/ml alone and in combination with rhSCF and rhIL-3 (R&D Systems, Minneapolis, MN, Cat# 203-ML) at 30 ng/ml. After one hour, 10  $\mu$ l of prepared cytokines, 50  $\mu$ l of the supernatants prepared in Example 30 (supernatants at 1:2 dilution = 50  $\mu$ l) and 20  $\mu$ l of diluted cells are added to the media which is already present in the wells to allow for a final total volume of 100  $\mu$ l. The plates are then placed in a 37°C/5% CO<sub>2</sub> incubator for five days.

[01214] Eighteen hours before the assay is harvested, 0.5  $\mu$ Ci/well of [3H] Thymidine is added in a 10  $\mu$ l volume to each well to determine the proliferation rate. The experiment is terminated by harvesting the cells from each 96-well plate to a filtermat

using the Tomtec Harvester 96. After harvesting, the filtermats are dried, trimmed and placed into OmniFilter assemblies consisting of one OmniFilter plate and one OmniFilter Tray. 60 µl Microscint is added to each well and the plate sealed with TopSeal-A press-on sealing film. A bar code 15 sticker is affixed to the first plate for counting. The sealed plates are then loaded and the level of radioactivity determined via the Packard Top Count and the printed data collected for analysis. The level of radioactivity reflects the amount of cell proliferation.

[01215] The studies described in this example test the activity of a given polypeptide to stimulate bone marrow CD34+ cell proliferation. One skilled in the art could easily modify the exemplified studies to test the activity of polynucleotides (e.g., gene therapy), antibodies, agonists, and/or antagonists and fragments and variants thereof. As a nonlimiting example, potential antagonists tested in this assay would be expected to inhibit cell proliferation in the presence of cytokines and/or to increase the inhibition of cell proliferation in the presence of cytokines and a given polypeptide. In contrast, potential agonists tested in this assay would be expected to enhance cell proliferation and/or to decrease the inhibition of cell proliferation in the presence of cytokines and a given polypeptide.

[01216] The ability of a gene to stimulate the proliferation of bone marrow CD34+ cells indicates that polynucleotides and polypeptides corresponding to the gene are useful for the diagnosis and treatment of disorders affecting the immune system and hematopoiesis. Representative uses are described in the "Immune Activity" and "Infectious Disease" sections above, and elsewhere herein.

***Example 41: Assay for Extracellular Matrix Enhanced Cell Response  
(EMECCR)***

[01217] The objective of the Extracellular Matrix Enhanced Cell Response (EMECCR) assay is to identify gene products (e.g., isolated polypeptides) that act on the hematopoietic stem cells in the context of the extracellular matrix (ECM) induced signal.

[01218] Cells respond to the regulatory factors in the context of signal(s) received from the surrounding microenvironment. For example, fibroblasts, and endothelial and epithelial stem cells fail to replicate in the absence of signals from the ECM.

Hematopoietic stem cells can undergo self-renewal in the bone marrow, but not in *in vitro* suspension culture. The ability of stem cells to undergo self-renewal *in vitro* is dependent upon their interaction with the stromal cells and the ECM protein fibronectin (fn). Adhesion of cells to fn is mediated by the  $\alpha_5\beta_1$  and  $\alpha_4\beta_1$  integrin receptors, which are expressed by human and mouse hematopoietic stem cells. The factor(s) which integrate with the ECM environment and are responsible for stimulating stem cell self-renewal have not yet been identified. Discovery of such factors should be of great interest in gene therapy and bone marrow transplant applications

[01219] Briefly, polystyrene, non tissue culture treated, 96-well plates are coated with fn fragment at a coating concentration of  $0.2 \mu\text{g}/\text{cm}^2$ . Mouse bone marrow cells are plated (1,000 cells/well) in 0.2 ml of serum-free medium. Cells cultured in the presence of IL-3 (5 ng/ml) + SCF (50 ng/ml) would serve as the positive control, conditions under which little self-renewal but pronounced differentiation of the stem cells is to be expected. Gene products of the invention (e.g., including, but not limited to, polynucleotides and polypeptides of the present invention, and supernatants produced in Example 30), are tested with appropriate negative controls in the presence and absence of SCF(5.0 ng/ml), where test factor supernatants represent 10% of the total assay volume. The plated cells are then allowed to grow by incubating in a low oxygen environment (5%  $\text{CO}_2$ , 7%  $\text{O}_2$ , and 88%  $\text{N}_2$ ) tissue culture incubator for 7 days. The number of proliferating cells within the wells is then quantitated by measuring thymidine incorporation into cellular DNA. Verification of the positive hits in the assay will require phenotypic characterization of the cells, which can be accomplished by scaling up of the culture system and using appropriate antibody reagents against cell surface antigens and FACScan.

[01220] One skilled in the art could easily modify the exemplified studies to test the activity of polynucleotides (e.g., gene therapy), antibodies, agonists, and/or antagonists and fragments and variants thereof.

[01221] If a particular polypeptide of the present invention is found to be a stimulator of hematopoietic progenitors, polynucleotides and polypeptides corresponding to the gene encoding said polypeptide may be useful for the diagnosis and treatment of disorders affecting the immune system and hematopoiesis. Representative uses are described in the "Immune Activity" and "Infectious Disease" sections above, and

elsewhere herein. The gene product may also be useful in the expansion of stem cells and committed progenitors of various blood lineages, and in the differentiation and/or proliferation of various cell types.

[01222] Additionally, the polynucleotides and/or polypeptides of the gene of interest and/or agonists and/or antagonists thereof, may also be employed to inhibit the proliferation and differentiation of hematopoietic cells and therefore may be employed to protect bone marrow stem cells from chemotherapeutic agents during chemotherapy. This antiproliferative effect may allow administration of higher doses of chemotherapeutic agents and, therefore, more effective chemotherapeutic treatment.

[01223] Moreover, polynucleotides and polypeptides corresponding to the gene of interest may also be useful for the treatment and diagnosis of hematopoietic related disorders such as, for example, anemia, pancytopenia, leukopenia, thrombocytopenia or leukemia since stromal cells are important in the production of cells of hematopoietic lineages. The uses include bone marrow cell ex-vivo culture, bone marrow transplantation, bone marrow reconstitution, radiotherapy or chemotherapy of neoplasia.

***Example 42: Human Dermal Fibroblast and Aortic Smooth Muscle Cell Proliferation***

[01224] The polypeptide of interest is added to cultures of normal human dermal fibroblasts (NHDF) and human aortic smooth muscle cells (AoSMC) and two co-assays are performed with each sample. The first assay examines the effect of the polypeptide of interest on the proliferation of normal human dermal fibroblasts (NHDF) or aortic smooth muscle cells (AoSMC). Aberrant growth of fibroblasts or smooth muscle cells is a part of several pathological processes, including fibrosis, and restenosis. The second assay examines IL6 production by both NHDF and SMC. IL6 production is an indication of functional activation. Activated cells will have increased production of a number of cytokines and other factors, which can result in a proinflammatory or immunomodulatory outcome. Assays are run with and without co-TNF $\alpha$  stimulation, in order to check for costimulatory or inhibitory activity.

[01225] Briefly, on day 1, 96-well black plates are set up with 1000 cells/well (NHDF) or 2000 cells/well (AoSMC) in 100  $\mu$ l culture media. NHDF culture media

contains: Clonetics FB basal media, 1mg/ml hFGF, 5mg/ml insulin, 50mg/ml gentamycin, 2%FBS, while AoSMC culture media contains Clonetics SM basal media, 0.5 µg/ml hEGF, 5mg/ml insulin, 1µg/ml hFGF, 50mg/ml gentamycin, 50 µg/ml Amphotericin B, 5%FBS. After incubation at 37°C for at least 4-5 hours culture media is aspirated and replaced with growth arrest media. Growth arrest media for NHDF contains fibroblast basal media, 50mg/ml gentamycin, 2% FBS, while growth arrest media for AoSMC contains SM basal media, 50mg/ml gentamycin, 50µg/ml Amphotericin B, 0.4% FBS. Incubate at 37 °C until day 2.

**[01226]** On day 2, serial dilutions and templates of the polypeptide of interest are designed such that they always include media controls and known-protein controls. For both stimulation and inhibition experiments, proteins are diluted in growth arrest media. For inhibition experiments, TNFa is added to a final concentration of 2ng/ml (NHDF) or 5ng/ml (AoSMC). Add 1/3 vol media containing controls or polypeptides of the present invention and incubate at 37 degrees C/5% CO<sub>2</sub> until day 5.

**[01227]** Transfer 60µl from each well to another labeled 96-well plate, cover with a plate-sealer, and store at 4 degrees C until Day 6 (for IL6 ELISA). To the remaining 100 µl in the cell culture plate, aseptically add Alamar Blue in an amount equal to 10% of the culture volume (10µl). Return plates to incubator for 3 to 4 hours. Then measure fluorescence with excitation at 530nm and emission at 590nm using the CytoFluor. This yields the growth stimulation/inhibition data.

**[01228]** On day 5, the IL6 ELISA is performed by coating a 96 well plate with 50-100 ul/well of Anti-Human IL6 Monoclonal antibody diluted in PBS, pH 7.4, incubate ON at room temperature.

**[01229]** On day 6, empty the plates into the sink and blot on paper towels. Prepare Assay Buffer containing PBS with 4% BSA. Block the plates with 200 µl/well of Pierce Super Block blocking buffer in PBS for 1-2 hr and then wash plates with wash buffer (PBS, 0.05% Tween-20). Blot plates on paper towels. Then add 50 µl/well of diluted Anti-Human IL-6 Monoclonal, Biotin-labeled antibody at 0.50 mg/ml. Make dilutions of IL-6 stock in media (30, 10, 3, 1, 0.3, 0 ng/ml). Add duplicate samples to top row of plate. Cover the plates and incubate for 2 hours at RT on shaker.

**[01230]** Plates are washed with wash buffer and blotted on paper towels. Dilute EU-labeled Streptavidin 1:1000 in Assay buffer, and add 100 µl/well. Cover the plate and

incubate 1 h at RT. Plates are again washed with wash buffer and blotted on paper towels.

[01231] Add 100  $\mu$ l/well of Enhancement Solution. Shake for 5 minutes. Read the plate on the Wallac DELFIA Fluorometer. Readings from triplicate samples in each assay were tabulated and averaged.

[01232] A positive result in this assay suggests AoSMC cell proliferation and that the polypeptide of the present invention may be involved in dermal fibroblast proliferation and/or smooth muscle cell proliferation. A positive result also suggests many potential uses of polypeptides, polynucleotides, agonists and/or antagonists of the polynucleotide/polypeptide of the present invention which gives a positive result. For example, inflammation and immune responses, wound healing, and angiogenesis, as detailed throughout this specification. Particularly, polypeptides of the present invention and polynucleotides of the present invention may be used in wound healing and dermal regeneration, as well as the promotion of vasculogenesis, both of the blood vessels and lymphatics. The growth of vessels can be used in the treatment of, for example, cardiovascular diseases. Additionally, antagonists of polypeptides and polynucleotides of the invention may be useful in treating diseases, disorders, and/or conditions which involve angiogenesis by acting as an anti-vascular agent (e.g., anti-angiogenesis). These diseases, disorders, and/or conditions are known in the art and/or are described herein, such as, for example, malignancies, solid tumors, benign tumors, for example hemangiomas, acoustic neuromas, neurofibromas, trachomas, and pyogenic granulomas; arteriosclerotic plaques; ocular angiogenic diseases, for example, diabetic retinopathy, retinopathy of prematurity, macular degeneration, corneal graft rejection, neovascular glaucoma, retrolental fibroplasia, rubeosis, retinoblastoma, uveitis and Pterygia (abnormal blood vessel growth) of the eye; rheumatoid arthritis; psoriasis; delayed wound healing; endometriosis; vasculogenesis; granulations; hypertrophic scars (keloids); nonunion fractures; scleroderma; trachoma; vascular adhesions; myocardial angiogenesis; coronary collaterals; cerebral collaterals; arteriovenous malformations; ischemic limb angiogenesis; Osler-Webber Syndrome; plaque neovascularization; telangiectasia; hemophilic joints; angiofibroma; fibromuscular dysplasia; wound granulation; Crohn's disease; and atherosclerosis. Moreover, antagonists of polypeptides



and polynucleotides of the invention may be useful in treating anti-hyperproliferative diseases and/or anti-inflammatory known in the art and/or described herein.

[01233] One skilled in the art could easily modify the exemplified studies to test the activity of polynucleotides (e.g., gene therapy), antibodies, agonists, and/or antagonists and fragments and variants thereof.

***Example 43: Cellular Adhesion Molecule (CAM) Expression  
on Endothelial Cells***

[01234] The recruitment of lymphocytes to areas of inflammation and angiogenesis involves specific receptor-ligand interactions between cell surface adhesion molecules (CAMs) on lymphocytes and the vascular endothelium. The adhesion process, in both normal and pathological settings, follows a multi-step cascade that involves intercellular adhesion molecule-1 (ICAM-1), vascular cell adhesion molecule-1 (VCAM-1), and endothelial leukocyte adhesion molecule-1 (E-selectin) expression on endothelial cells (EC). The expression of these molecules and others on the vascular endothelium determines the efficiency with which leukocytes may adhere to the local vasculature and extravasate into the local tissue during the development of an inflammatory response. The local concentration of cytokines and growth factor participate in the modulation of the expression of these CAMs.

[01235] Briefly, endothelial cells (e.g., Human Umbilical Vein Endothelial cells (HUVECs)) are grown in a standard 96 well plate to confluence, growth medium is removed from the cells and replaced with 100  $\mu$ l of 199 Medium (10% fetal bovine serum (FBS)). Samples for testing and positive or negative controls are added to the plate in triplicate (in 10  $\mu$ l volumes). Plates are then incubated at 37°C for either 5 h (selectin and integrin expression) or 24 h (integrin expression only). Plates are aspirated to remove medium and 100  $\mu$ l of 0.1% paraformaldehyde-PBS(with Ca<sup>++</sup> and Mg<sup>++</sup>) is added to each well. Plates are held at 4°C for 30 min. Fixative is removed from the wells and wells are washed 1X with PBS(+Ca,Mg) + 0.5% BSA and drained. 10  $\mu$ l of diluted primary antibody is added to the test and control wells. Anti-ICAM-1-Biotin, Anti-VCAM-1-Biotin and Anti-E-selectin-Biotin are used at a concentration of 10  $\mu$ g/ml

(1:10 dilution of 0.1 mg/ml stock antibody). Cells are incubated at 37°C for 30 min. in a humidified environment. Wells are washed three times with PBS(+Ca,Mg) + 0.5% BSA. 20 µl of diluted ExtrAvidin-Alkaline Phosphatase (1:5,000 dilution, referred to herein as the working dilution) are added to each well and incubated at 37°C for 30 min. Wells are washed three times with PBS(+Ca,Mg)+0.5% BSA. Dissolve 1 tablet of p-Nitrophenol Phosphate pNPP per 5 ml of glycine buffer (pH 10.4). 100 µl of pNPP substrate in glycine buffer is added to each test well. Standard wells in triplicate are prepared from the working dilution of the ExtrAvidin-Alkaline Phosphatase in glycine buffer: 1:5,000 ( $10^0$ ) >  $10^{-0.5}$  >  $10^{-1}$  >  $10^{-1.5}$ . 5 µl of each dilution is added to triplicate wells and the resulting AP content in each well is 5.50 ng, 1.74 ng, 0.55 ng, 0.18 ng. 100 µl of pNPP reagent is then added to each of the standard wells. The plate is incubated at 37°C for 4h. A volume of 50 µl of 3M NaOH is added to all wells. The plate is read on a plate reader at 405 nm using the background subtraction option on blank wells filled with glycine buffer only. Additionally, the template is set up to indicate the concentration of AP-conjugate in each standard well [5.50 ng; 1.74 ng; 0.55 ng; 0.18 ng]. Results are indicated as amount of bound AP-conjugate in each sample.

#### ***Example 44: Alamar Blue Endothelial Cells Proliferation Assay***

[01236] This assay may be used to quantitatively determine protein mediated inhibition of bFGF-induced proliferation of Bovine Lymphatic Endothelial Cells (LECs), Bovine Aortic Endothelial Cells (BAECs) or Human Microvascular Uterine Myometrial Cells (UTMECs). This assay incorporates a fluorometric growth indicator based on detection of metabolic activity. A standard Alamar Blue Proliferation Assay is prepared in EGM-2MV with 10 ng /ml of bFGF added as a source of endothelial cell stimulation. This assay may be used with a variety of endothelial cells with slight changes in growth medium and cell concentration. Dilutions of the protein batches to be tested are diluted as appropriate. Serum-free medium (GIBCO SFM) without bFGF is used as a non-stimulated control and Angiostatin or TSP-1 are included as a known inhibitory controls.

[01237] Briefly, LEC, BAECs or UTMECs are seeded in growth media at a density of 5000 to 2000 cells/well in a 96 well plate and placed at 37degrees C overnight. After the overnight incubation of the cells, the growth media is removed and replaced with

GIBCO EC-SFM. The cells are treated with the appropriate dilutions of the protein of interest or control protein sample(s) (prepared in SFM ) in triplicate wells with additional bFGF to a concentration of 10 ng/ ml. Once the cells have been treated with the samples, the plate(s) is/are placed back in the 37° C incubator for three days. After three days 10 ml of stock alamar blue (Biosource Cat# DAL1100) is added to each well and the plate(s) is/are placed back in the 37°C incubator for four hours. The plate(s) are then read at 530nm excitation and 590nm emission using the CytoFluor fluorescence reader. Direct output is recorded in relative fluorescence units.

[01238] Alamar blue is an oxidation-reduction indicator that both fluoresces and changes color in response to chemical reduction of growth medium resulting from cell growth. As cells grow in culture, innate metabolic activity results in a chemical reduction of the immediate surrounding environment. Reduction related to growth causes the indicator to change from oxidized (non-fluorescent blue) form to reduced (fluorescent red) form (i.e., stimulated proliferation will produce a stronger signal and inhibited proliferation will produce a weaker signal and the total signal is proportional to the total number of cells as well as their metabolic activity). The background level of activity is observed with the starvation medium alone. This is compared to the output observed from the positive control samples (bFGF in growth medium) and protein dilutions.

#### ***Example 45: Detection of Inhibition of a Mixed Lymphocyte Reaction***

[01239] This assay can be used to detect and evaluate inhibition of a Mixed Lymphocyte Reaction (MLR) by gene products (e.g., isolated polypeptides). Inhibition of a MLR may be due to a direct effect on cell proliferation and viability, modulation of costimulatory molecules on interacting cells, modulation of adhesiveness between lymphocytes and accessory cells, or modulation of cytokine production by accessory cells. Multiple cells may be targeted by these polypeptides since the peripheral blood mononuclear fraction used in this assay includes T, B and natural killer lymphocytes, as well as monocytes and dendritic cells.

[01240] Polypeptides of interest found to inhibit the MLR may find application in diseases associated with lymphocyte and monocyte activation or proliferation. These include, but are not limited to, diseases such as asthma, arthritis, diabetes, inflammatory

skin conditions, psoriasis, eczema, systemic lupus erythematosus, multiple sclerosis, glomerulonephritis, inflammatory bowel disease, crohn's disease, ulcerative colitis, arteriosclerosis, cirrhosis, graft vs. host disease, host vs. graft disease, hepatitis, leukemia and lymphoma.

[01241] Briefly, PBMCs from human donors are purified by density gradient centrifugation using Lymphocyte Separation Medium (LSM<sup>®</sup>, density 1.0770 g/ml, Organon Teknika Corporation, West Chester, PA). PBMCs from two donors are adjusted to  $2 \times 10^6$  cells/ml in RPMI-1640 (Life Technologies, Grand Island, NY) supplemented with 10% FCS and 2 mM glutamine. PBMCs from a third donor is adjusted to  $2 \times 10^5$  cells/ml. Fifty microliters of PBMCs from each donor is added to wells of a 96-well round bottom microtiter plate. Dilutions of test materials (50  $\mu$ l) is added in triplicate to microtiter wells. Test samples (of the protein of interest) are added for final dilution of 1:4; rhuIL-2 (R&D Systems, Minneapolis, MN, catalog number 202-IL) is added to a final concentration of 1  $\mu$ g/ml; anti-CD4 mAb (R&D Systems, clone 34930.11, catalog number MAB379) is added to a final concentration of 10  $\mu$ g/ml. Cells are cultured for 7-8 days at 37°C in 5% CO<sub>2</sub>, and 1  $\mu$ C of [<sup>3</sup>H] thymidine is added to wells for the last 16 hrs of culture. Cells are harvested and thymidine incorporation determined using a Packard TopCount. Data is expressed as the mean and standard deviation of triplicate determinations.

[01242] Samples of the protein of interest are screened in separate experiments and compared to the negative control treatment, anti-CD4 mAb, which inhibits proliferation of lymphocytes and the positive control treatment, IL-2 (either as recombinant material or supernatant), which enhances proliferation of lymphocytes.

[01243] One skilled in the art could easily modify the exemplified studies to test the activity of polynucleotides (e.g., gene therapy), antibodies, agonists, and/or antagonists and fragments and variants thereof.

#### ***Example 46: Assays for Protease Activity***

[01244] The following assay may be used to assess protease activity of the polypeptides of the invention.

[01245] Gelatin and casein zymography are performed essentially as described (Heusen et al., Anal. Biochem., 102:196-202 (1980); Wilson et al., Journal of Urology, 149:653-658 (1993)). Samples are run on 10% polyacrylamide/0.1% SDS gels containing 1% gelatin or casein, soaked in 2.5% triton at room temperature for 1 hour, and in 0.1M glycine, pH 8.3 at 37°C 5 to 16 hours. After staining in amido black areas of proteolysis appear as clear areas against the blue-black background. Trypsin (Sigma T8642) is used as a positive control.

[01246] Protease activity is also determined by monitoring the cleavage of n-a-benzoyl-L-arginine ethyl ester (BAEE) (Sigma B-4500). Reactions are set up in (25mM NaPO<sub>4</sub>, 1mM EDTA, and 1mM BAEE), pH 7.5. Samples are added and the change in absorbance at 260nm is monitored on the Beckman DU-6 spectrophotometer in the time-drive mode. Trypsin is used as a positive control.

[01247] Additional assays based upon the release of acid-soluble peptides from casein or hemoglobin measured as absorbance at 280 nm or colorimetrically using the Folin method are performed as described in Bergmeyer, et al., Methods of Enzymatic Analysis, 5 (1984). Other assays involve the solubilization of chromogenic substrates (Ward, Applied Science, 251-317 (1983)).

#### ***Example 47: Identifying Serine Protease Substrate Specificity***

[01248] Methods known in the art or described herein may be used to determine the substrate specificity of the polypeptides of the present invention having serine protease activity. A preferred method of determining substrate specificity is by the use of positional scanning synthetic combinatorial libraries as described in GB 2 324 529 (incorporated herein in its entirety).

#### ***Example 48: Ligand Binding Assays***

[01249] The following assay may be used to assess ligand binding activity of the polypeptides of the invention.

[01250] Ligand binding assays provide a direct method for ascertaining receptor pharmacology and are adaptable to a high throughput format. The purified ligand for a

polypeptide is radiolabeled to high specific activity (50-2000 Ci/mmol) for binding studies. A determination is then made that the process of radiolabeling does not diminish the activity of the ligand towards its polypeptide. Assay conditions for buffers, ions, pH and other modulators such as nucleotides are optimized to establish a workable signal to noise ratio for both membrane and whole cell polypeptide sources. For these assays, specific polypeptide binding is defined as total associated radioactivity minus the radioactivity measured in the presence of an excess of unlabeled competing ligand. Where possible, more than one competing ligand is used to define residual nonspecific binding.

#### ***Example 49: Functional Assay in *Xenopus* Oocytes***

[01251] Capped RNA transcripts from linearized plasmid templates encoding the polypeptides of the invention are synthesized *in vitro* with RNA polymerases in accordance with standard procedures. *In vitro* transcripts are suspended in water at a final concentration of 0.2 mg/ml. Ovarian lobes are removed from adult female toads, Stage V defolliculated oocytes are obtained, and RNA transcripts (10 ng/oocyte) are injected in a 50 nl bolus using a microinjection apparatus. Two electrode voltage clamps are used to measure the currents from individual *Xenopus oocytes* in response polypeptides and polypeptide agonist exposure. Recordings are made in Ca<sup>2+</sup> free Barth's medium at room temperature. The *Xenopus* system can be used to screen known ligands and tissue/cell extracts for activating ligands.

#### ***Example 50: Microphysiometric Assays***

[01252] Activation of a wide variety of secondary messenger systems results in extrusion of small amounts of acid from a cell. The acid formed is largely as a result of the increased metabolic activity required to fuel the intracellular signaling process. The pH changes in the media surrounding the cell are very small but are detectable by the CYTOSENSOR microphysiometer (Molecular Devices Ltd., Menlo Park, Calif.). The CYTOSENSOR is thus capable of detecting the activation of polypeptide which is coupled to an energy utilizing intracellular signaling pathway.

***Example 51: Extract/Cell Supernatant Screening***

[01253] A large number of mammalian receptors exist for which there remains, as yet, no cognate activating ligand (agonist). Thus, active ligands for these receptors may not be included within the ligands banks as identified to date. Accordingly, the polypeptides of the invention can also be functionally screened (using calcium, cAMP, microphysiometer, oocyte electrophysiology, etc., functional screens) against tissue extracts to identify its natural ligands. Extracts that produce positive functional responses can be sequentially subfractionated until an activating ligand is isolated and identified.

***Example 52: Calcium and cAMP Functional Assays***

[01254] Seven transmembrane receptors which are expressed in HEK 293 cells have been shown to be coupled functionally to activation of PLC and calcium mobilization and/or cAMP stimulation or inhibition. Basal calcium levels in the HEK 293 cells in receptor-transfected or vector control cells were observed to be in the normal, 100 nM to 200 nM, range. HEK 293 cells expressing recombinant receptors are loaded with fura 2 and in a single day >150 selected ligands or tissue/cell extracts are evaluated for agonist induced calcium mobilization. Similarly, HEK 293 cells expressing recombinant receptors are evaluated for the stimulation or inhibition of cAMP production using standard cAMP quantitation assays. Agonists presenting a calcium transient or cAMP fluctuation are tested in vector control cells to determine if the response is unique to the transfected cells expressing receptor.

***Example 53: ATP-binding assay***

[01255] The following assay may be used to assess ATP-binding activity of polypeptides of the invention.

[01256] ATP-binding activity of the polypeptides of the invention may be detected using the ATP-binding assay described in U.S. Patent 5,858,719, which is herein

incorporated by reference in its entirety. Briefly, ATP-binding to polypeptides of the invention is measured via photoaffinity labeling with 8-azido-ATP in a competition assay. Reaction mixtures containing 1 mg/ml of the ABC transport protein of the present invention are incubated with varying concentrations of ATP, or the non-hydrolyzable ATP analog adenylyl-5'-imidodiphosphate for 10 minutes at 4°C. A mixture of 8-azido-ATP (Sigma Chem. Corp., St. Louis, MO.) plus 8-azido-ATP ( $^{32}\text{P}$ -ATP) (5 mCi/ $\mu\text{mol}$ , ICN, Irvine CA.) is added to a final concentration of 100  $\mu\text{M}$  and 0.5 ml aliquots are placed in the wells of a porcelain spot plate on ice. The plate is irradiated using a short wave 254 nm UV lamp at a distance of 2.5 cm from the plate for two one-minute intervals with a one-minute cooling interval in between. The reaction is stopped by addition of dithiothreitol to a final concentration of 2mM. The incubations are subjected to SDS-PAGE electrophoresis, dried, and autoradiographed. Protein bands corresponding to the particular polypeptides of the invention are excised, and the radioactivity quantified. A decrease in radioactivity with increasing ATP or adenylyl-5'-imidodiphosphate provides a measure of ATP affinity to the polypeptides.

#### ***Example 54: Small Molecule Screening***

[01257] This invention is particularly useful for screening therapeutic compounds by using the polypeptides of the invention, or binding fragments thereof, in any of a variety of drug screening techniques. The polypeptide or fragment employed in such a test may be affixed to a solid support, expressed on a cell surface, free in solution, or located intracellularly. One method of drug screening utilizes eukaryotic or prokaryotic host cells which are stably transformed with recombinant nucleic acids expressing the polypeptide or fragment. Drugs are screened against such transformed cells in competitive binding assays. One may measure, for example, the formulation of complexes between the agent being tested and polypeptide of the invention.

[01258] Thus, the present invention provides methods of screening for drugs or any other agents which affect activities mediated by the polypeptides of the invention. These methods comprise contacting such an agent with a polypeptide of the invention or fragment thereof and assaying for the presence of a complex between the agent and the



polypeptide or fragment thereof, by methods well known in the art. In such a competitive binding assay, the agents to screen are typically labeled. Following incubation, free agent is separated from that present in bound form, and the amount of free or uncomplexed label is a measure of the ability of a particular agent to bind to the polypeptides of the invention.

[01259] Another technique for drug screening provides high throughput screening for compounds having suitable binding affinity to the polypeptides of the invention, and is described in great detail in European Patent Application 84/03564, published on September 13, 1984, which is herein incorporated by reference in its entirety. Briefly stated, large numbers of different small molecule test compounds are synthesized on a solid substrate, such as plastic pins or some other surface. The test compounds are reacted with polypeptides of the invention and washed. Bound polypeptides are then detected by methods well known in the art. Purified polypeptides are coated directly onto plates for use in the aforementioned drug screening techniques. In addition, non-neutralizing antibodies may be used to capture the peptide and immobilize it on the solid support.

[01260] This invention also contemplates the use of competitive drug screening assays in which neutralizing antibodies capable of binding polypeptides of the invention specifically compete with a test compound for binding to the polypeptides or fragments thereof. In this manner, the antibodies are used to detect the presence of any peptide which shares one or more antigenic epitopes with a polypeptide of the invention.

### ***Example 55: Phosphorylation Assay***

[01261] In order to assay for phosphorylation activity of the polypeptides of the invention, a phosphorylation assay as described in U.S. Patent 5,958,405 (which is herein incorporated by reference) is utilized. Briefly, phosphorylation activity may be measured by phosphorylation of a protein substrate using gamma-labeled  $^{32}\text{P}$ -ATP and quantitation of the incorporated radioactivity using a gamma radioisotope counter. The polypeptides of the invention are incubated with the protein substrate,  $^{32}\text{P}$ -ATP, and a kinase buffer. The  $^{32}\text{P}$  incorporated into the substrate is then separated from free  $^{32}\text{P}$ -ATP by electrophoresis, and the incorporated  $^{32}\text{P}$  is counted and compared to a negative control.

Radioactivity counts above the negative control are indicative of phosphorylation activity of the polypeptides of the invention.

***Example 56: Detection of Phosphorylation Activity (Activation) of the Polypeptides of the Invention in the Presence of Polypeptide Ligands***

[01262] Methods known in the art or described herein may be used to determine the phosphorylation activity of the polypeptides of the invention. A preferred method of determining phosphorylation activity is by the use of the tyrosine phosphorylation assay as described in U.S. Patent 5,817,471 (incorporated herein by reference).

***Example 57: Identification Of Signal Transduction Proteins That Interact With Polypeptides Of The Present Invention***

[01263] The purified polypeptides of the invention are research tools for the identification, characterization and purification of additional signal transduction pathway proteins or receptor proteins. Briefly, labeled receptor PTK polypeptide is useful as a reagent for the purification of molecules with which it interacts. In one embodiment of affinity purification, receptor PTK polypeptide is covalently coupled to a chromatography column. Cell-free extract derived from putative target cells, such as carcinoma tissues, is passed over the column, and molecules with appropriate affinity bind to the receptor PTK polypeptides, or specific phosphotyrosine-recognition domains thereof. The receptor PTK polypeptide interacting protein-complex is recovered from the column, dissociated, and the recovered molecule subjected to N-terminal protein sequencing. This amino acid sequence is then used to identify the captured molecule or to design degenerate oligonucleotide probes for cloning the relevant gene from an appropriate cDNA library.

***Example 58: IL-6 Bioassay***

[01264] To test the proliferative effects of the polypeptides of the invention, the IL-6 Bioassay as described by Marz *et al.* is utilized (*Proc. Natl. Acad. Sci., U.S.A.*, 95:3251-56 (1998), which is herein incorporated by reference). Briefly, IL-6 dependent B9 murine cells are washed three times in IL-6 free medium and plated at a concentration of 5,000 cells per well in 50  $\mu$ l, and 50  $\mu$ l of the IL-6-like polypeptide is added. After 68 hrs. at 37°C, the number of viable cells is measured by adding the tetrazolium salt thiazolyl blue (MTT) and incubating for a further 4 hrs. at 37°C. B9 cells are lysed by SDS and optical density is measured at 570 nm. Controls containing IL-6 (positive) and no cytokine (negative) are utilized. Enhanced proliferation in the test sample(s) relative to the negative control is indicative of proliferative effects mediated by polypeptides of the invention.

#### ***Example 59: Support of Chicken Embryo Neuron Survival***

[01265] To test whether sympathetic neuronal cell viability is supported by polypeptides of the invention, the chicken embryo neuronal survival assay of Senaldi *et al* is utilized (*Proc. Natl. Acad. Sci., U.S.A.*, 96:11458-63 (1998), which is herein incorporated by reference). Briefly, motor and sympathetic neurons are isolated from chicken embryos, resuspended in L15 medium (with 10% FCS, glucose, sodium selenite, progesterone, conalbumin, putrescine, and insulin; Life Technologies, Rockville, MD.) and Dulbecco's modified Eagles medium [with 10% FCS, glutamine, penicillin, and 25 mM Hepes buffer (pH 7.2); Life Technologies, Rockville, MD.], respectively, and incubated at 37°C in 5% CO<sub>2</sub> in the presence of different concentrations of the purified IL-6-like polypeptide, as well as a negative control lacking any cytokine. After 3 days, neuron survival is determined by evaluation of cellular morphology, and through the use of the colorimetric assay of Mosmann (Mossmann, T., J. Immunol. Methods, 65:55-63 (1983)). Enhanced neuronal cell viability as compared to the controls lacking cytokine is indicative of the ability of the inventive purified IL-6-like polypeptide(s) to enhance the survival of neuronal cells.

#### ***Example 60: Assay for Phosphatase Activity***

[01266] The following assay may be used to assess serine/threonine phosphatase (PTPase) activity of the polypeptides of the invention.

[01267] In order to assay for serine/threonine phosphatase (PTPase) activity, assays can be utilized which are widely known to those skilled in the art. For example, the serine/threonine phosphatase (PSPase) activity is measured using a PSPase assay kit from New England Biolabs, Inc. Myelin basic protein (MyBP), a substrate for PSPase, is phosphorylated on serine and threonine residues with cAMP-dependent Protein Kinase in the presence of [ $^{32}$ P]ATP. Protein serine/threonine phosphatase activity is then determined by measuring the release of inorganic phosphate from  $^{32}$ P-labeled MyBP.

***Example 61: Interaction of Serine/Threonine Phosphatases with other Proteins***

[01268] The polypeptides of the invention with serine/threonine phosphatase activity as determined in Example 60 are research tools for the identification, characterization and purification of additional interacting proteins or receptor proteins, or other signal transduction pathway proteins. Briefly, labeled polypeptide(s) of the invention is useful as a reagent for the purification of molecules with which it interacts. In one embodiment of affinity purification, polypeptide of the invention is covalently coupled to a chromatography column. Cell-free extract derived from putative target cells, such as neural or liver cells, is passed over the column, and molecules with appropriate affinity bind to the polypeptides of the invention. The polypeptides of the invention -complex is recovered from the column, dissociated, and the recovered molecule subjected to N-terminal protein sequencing. This amino acid sequence is then used to identify the captured molecule or to design degenerate oligonucleotide probes for cloning the relevant gene from an appropriate cDNA library.

***Example 62: Assaying for Heparanase Activity***

[01269] In order to assay for heparanase activity of the polypeptides of the invention, the heparanase assay described by Vlodavsky et al is utilized (Vlodavsky, I., et al., Nat. Med., 5:793-802 (1999)). Briefly, cell lysates, conditioned media or intact cells ( $1 \times 10^6$

cells per 35-mm dish) are incubated for 18 hrs at 37°C, pH 6.2-6.6, with <sup>35</sup>S-labeled ECM or soluble ECM derived peak I proteoglycans. The incubation medium is centrifuged and the supernatant is analyzed by gel filtration on a Sepharose CL-6B column (0.9 x 30 cm). Fractions are eluted with PBS and their radioactivity is measured. Degradation fragments of heparan sulfate side chains are eluted from Sepharose 6B at  $0.5 < K_{av} < 0.8$  (peak II). Each experiment is done at least three times. Degradation fragments corresponding to "peak II," as described by Vlodavsky et al., is indicative of the activity of the polypeptides of the invention in cleaving heparan sulfate.

### ***Example 63: Immobilization of biomolecules***

[01270] This example provides a method for the stabilization of polypeptides of the invention in non-host cell lipid bilayer constructs (see, e.g., Bieri et al., Nature Biotech 17:1105-1108 (1999), hereby incorporated by reference in its entirety herein) which can be adapted for the study of polypeptides of the invention in the various functional assays described above. Briefly, carbohydrate-specific chemistry for biotinylation is used to confine a biotin tag to the extracellular domain of the polypeptides of the invention, thus allowing uniform orientation upon immobilization. A 50uM solution of polypeptides of the invention in washed membranes is incubated with 20 mM NaIO<sub>4</sub> and 1.5 mg/ml (4mM) BACH or 2 mg/ml (7.5mM) biotin-hydrazide for 1 hr at room temperature (reaction volume, 150ul). Then the sample is dialyzed (Pierce Slidealizer Cassett, 10 kDa cutoff; Pierce Chemical Co., Rockford IL) at 4C first for 5 h, exchanging the buffer after each hour, and finally for 12 h against 500 ml buffer R (0.15 M NaCl, 1 mM MgCl<sub>2</sub>, 10 mM sodium phosphate, pH7). Just before addition into a cuvette, the sample is diluted 1:5 in buffer ROG50 (Buffer R supplemented with 50 mM octylglucoside).

### ***Example 64: TAQMAN***

[01271] Quantitative PCR (QPCR). Total RNA from cells in culture are extracted by Trizol separation as recommended by the supplier (LifeTechnologies). (Total RNA is

treated with DNase I (Life Technologies) to remove any contaminating genomic DNA before reverse transcription.) Total RNA (50 ng) is used in a one-step, 50ul, RT-QPCR, consisting of Taqman Buffer A (Perkin-Elmer; 50 mM KCl/10 mM Tris, pH 8.3), 5.5 mM MgCl<sub>2</sub>, 240 μM each dNTP, 0.4 units RNase inhibitor(Promega), 8%glycerol, 0.012% Tween-20, 0.05% gelatin, 0.3uM primers, 0.1uM probe, 0.025units Amplitaq Gold (Perkin-Elmer) and 2.5 units Superscript II reverse transcriptase (Life Technologies). As a control for genomic contamination, parallel reactions are setup without reverse transcriptase. The relative abundance of (unknown) and 18S RNAs are assessed by using the Applied Biosystems Prism 7700 Sequence Detection System (Livak, K. J., Flood, S. J., Marmaro, J., Giusti, W. & Deetz, K. (1995) PCR Methods Appl. 4, 357-362). Reactions are carried out at 48°C for 30 min, 95°C for 10 min, followed by 40 cycles of 95°C for 15s, 60°C for 1 min. Reactions are performed in triplicate.

[01272] Primers (f & r) and FRET probes sets are designed using Primer Express Software (Perkin-Elmer). Probes are labeled at the 5'-end with the reporter dye 6-FAM and on the 3'-end with the quencher dye TAMRA (Biosource International, Camarillo, CA or Perkin-Elmer).

### ***Example 65: Assays for Metalloproteinase Activity***

[01273] Metalloproteinases (EC 3.4.24.-) are peptide hydrolases which use metal ions, such as Zn<sup>2+</sup>, as the catalytic mechanism. Metalloproteinase activity of polypeptides of the present invention can be assayed according to the following methods.

#### ***Proteolysis of alpha-2-macroglobulin***

[01274] To confirm protease activity, purified polypeptides of the invention are mixed with the substrate alpha-2-macroglobulin (0.2 unit/ml; Boehringer Mannheim, Germany) in 1x assay buffer (50 mM HEPES, pH 7.5, 0.2 M NaCl, 10 mM CaCl<sub>2</sub>, 25 μM ZnCl<sub>2</sub> and 0.05% Brij-35) and incubated at 37°C for 1-5 days. Trypsin is used as positive control. Negative controls contain only alpha-2-macroglobulin in assay buffer. The samples are collected and boiled in SDS-PAGE sample buffer containing 5% 2-mercaptoethanol for 5-min, then loaded onto 8% SDS-polyacrylamide gel. After

electrophoresis the proteins are visualized by silver staining. Proteolysis is evident by the appearance of lower molecular weight bands as compared to the negative control.

***Inhibition of alpha-2-macroglobulin proteolysis by inhibitors of metalloproteinases***

[01275] Known metalloproteinase inhibitors (metal chelators (EDTA, EGTA, AND  $\text{HgCl}_2$ ), peptide metalloproteinase inhibitors (TIMP-1 and TIMP-2), and commercial small molecule MMP inhibitors) are used to characterize the proteolytic activity of polypeptides of the invention. The three synthetic MMP inhibitors used are: MMP inhibitor I, [ $\text{IC}_{50}$  = 1.0  $\mu\text{M}$  against MMP-1 and MMP-8;  $\text{IC}_{50}$  = 30  $\mu\text{M}$  against MMP-9;  $\text{IC}_{50}$  = 150  $\mu\text{M}$  against MMP-3]; MMP-3 (stromelysin-1) inhibitor I [ $\text{IC}_{50}$  = 5  $\mu\text{M}$  against MMP-3], and MMP-3 inhibitor II [ $\text{K}_i$  = 130 nM against MMP-3]; inhibitors available through Calbiochem, catalog # 444250, 444218, and 444225, respectively). Briefly, different concentrations of the small molecule MMP inhibitors are mixed with purified polypeptides of the invention (50 $\mu\text{g}/\text{ml}$ ) in 22.9  $\mu\text{l}$  of 1x HEPES buffer (50 mM HEPES, pH 7.5, 0.2 M NaCl, 10 mM  $\text{CaCl}_2$ , 25  $\mu\text{M}$   $\text{ZnCl}_2$  and 0.05%Brij-35) and incubated at room temperature (24 °C) for 2-hr, then 7.1  $\mu\text{l}$  of substrate alpha-2-macroglobulin (0.2 unit/ml) is added and incubated at 37°C for 20-hr. The reactions are stopped by adding 4x sample buffer and boiled immediately for 5 minutes. After SDS-PAGE, the protein bands are visualized by silver stain.

***Synthetic Fluorogenic Peptide Substrates Cleavage Assay***

[01276] The substrate specificity for polypeptides of the invention with demonstrated metalloproteinase activity can be determined using synthetic fluorogenic peptide substrates (purchased from BACHEM Bioscience Inc). Test substrates include, M-1985, M-2225, M-2105, M-2110, and M-2255. The first four are MMP substrates and the last one is a substrate of tumor necrosis factor- $\alpha$  (TNF- $\alpha$ ) converting enzyme (TACE). All the substrates are prepared in 1:1 dimethyl sulfoxide (DMSO) and water. The stock solutions are 50-500  $\mu\text{M}$ . Fluorescent assays are performed by using a Perkin Elmer LS 50B luminescence spectrometer equipped with a constant temperature water bath. The excitation  $\lambda$  is 328 nm and the emission  $\lambda$  is 393 nm. Briefly, the assay is carried out by incubating 176  $\mu\text{l}$  1x HEPES buffer (0.2 M NaCl, 10 mM  $\text{CaCl}_2$ , 0.05% Brij-35 and 50 mM HEPES, pH 7.5) with 4  $\mu\text{l}$  of substrate solution (50  $\mu\text{M}$ ) at 25 °C for 15 minutes,

and then adding 20  $\mu$ l of a purified polypeptide of the invention into the assay cuvet. The final concentration of substrate is 1  $\mu$ M. Initial hydrolysis rates are monitored for 30-min.

***Example 66: Characterization of the cDNA contained in a deposited plasmid***

[01277] The size of the cDNA insert contained in a deposited plasmid may be routinely determined using techniques known in the art, such as PCR amplification using synthetic primers hybridizable to the 3' and 5' ends of the cDNA sequence. For example, two primers of 17-30 nucleotides derived from each end of the cDNA (i.e., hybridizable to the absolute 5' nucleotide or the 3' nucleotide end of the sequence of SEQ ID NO:X, respectively) are synthesized and used to amplify the cDNA using the deposited cDNA plasmid as a template. The polymerase chain reaction is carried out under routine conditions, for instance, in 25  $\mu$ l of reaction mixture with 0.5  $\mu$ g of the above cDNA template. A convenient reaction mixture is 1.5-5 mM  $MgCl_2$ , 0.01% (w/v) gelatin, 20  $\mu$ M each of dATP, dCTP, dGTP, dTTP, 25 pmol of each primer and 0.25 Unit of Taq polymerase. Thirty five cycles of PCR (denaturation at 94 degree C for 1 min; annealing at 55 degree C for 1 min; elongation at 72 degree C for 1 min) are performed with a Perkin-Elmer Cetus automated thermal cycler. The amplified product is analyzed by agarose gel electrophoresis. The PCR product is verified to be the selected sequence by subcloning and sequencing the DNA product.

[01278] Use of the above methodologies and/or other methodologies known in the art generates fragments from the clone corresponding to the approximate fragments described in Table 8, below. Accordingly, Table 8 provides a physical characterization of certain clones encompassed by the invention. The first column provides the unique clone identifier, "Clone ID NO:Z", for cDNA clones of the invention, as described in Table 1A. The second column provides the approximate size of the cDNA insert contained in the corresponding cDNA clone.



**TABLE 8**

Clone ID NO:Z	cDNA Insert Size:
-----	-----
HCUHF41	500
-----	-----
HCUHS19	300
-----	-----
HCWBU49	400
-----	-----
HDPBS54	2100
-----	-----
HDPJI05	2900
-----	-----
HHEKP61	1500
-----	-----
HMAKA11	1700
-----	-----
HMSG27	1600
-----	-----
HMSJW12	3100
-----	-----
HMVBO34	1200
-----	-----
HMVDU16	1400
-----	-----
HNGBC50	1700
-----	-----
HNGEP64	1200
-----	-----
HNGEU35	700
-----	-----
HNGFK28	900
-----	-----
HNHGQ79	1000
-----	-----
HNHIA73	300
-----	-----
HTOAF15	1400
-----	-----
HTXJN04	1500
-----	-----
HTXMY23	1100
-----	-----
HWADF26	900
-----	-----
HYABP53	1500
-----	-----

HWBFY46	1800
-----	
HWBAJ13	1800
-----	
HTWFJ71	600
-----	
HTOJL36	1700
-----	
HTOHH78	2100
-----	
HTOFO18	2200
-----	
HTGAD16	1000
-----	
HSAWC55	700
-----	
HSAUQ81	400
-----	
HSATC79	700
-----	
HNHOF94	800
-----	
HNHFF79	300
-----	
HNHEQ03	700
-----	
HNHDY14	800
-----	
HNHDV95	700
-----	
HNHDD52	300
-----	
HNHCA27	500
-----	
HNHBI48	600
-----	
HNHBF47	1000
-----	
HNHAB36	900
-----	
HNGPN54	500
-----	
HNGPL75	1100
-----	
HNGMW65	400
-----	
HNGKC69	800
-----	
HNGJE23	1100
-----	

HNGFX75	1100
HNGFH30	2800
HNGCE81	1700
HNGCE40	800
HNGCC02	500
HNFCV02	900
HNEDF15	1300
HMWIQ64	1600
HMWIM05	600
HMWHB71	1200
HMWGQ48	700
HMSHD80	1100
HMSFR23	2000
HMOAD90	400
HMOAD60	200
HMMBT03	700
HMMAB96	800
HMAIH15	500
HMAIA55	1100
HLYCQ43	800
HLMMK61	500
HLMBX91	900
HLMBQ48	400
HJABQ75	1100
HHEPT05	2500

HHEPT03	1100
-----	
HDTBD67	2200
-----	
HDTAT80	2400
-----	
HDTAA31	1500
-----	
HDPFZ70	800
-----	
HCWHF76	700
-----	
HCUHH95	2200
-----	
HCUDQ81	300
-----	
HCFON45	1000
-----	
HBMBX72	900
-----	

[01279] It will be clear that the invention may be practiced otherwise than as particularly described in the foregoing description and examples. Numerous modifications and variations of the present invention are possible in light of the above teachings and, therefore, are within the scope of the appended claims.

[01280] The entire disclosure of each document cited (including patents, patent applications, journal articles, abstracts, laboratory manuals, books, or other disclosures) in the Background of the Invention, Detailed Description, and Examples is hereby incorporated herein by reference. In addition, the CD-R copy of the sequence listing submitted herewith and the corresponding computer readable form are both incorporated herein by reference in their entireties. The specification and Sequence Listing of each of the following U.S. applications are herein incorporated by reference in their entirety: Application No. 60/179,065, filed on 31-Jan-2000; Application No. 60/180,628, filed on 04-Feb-2000; Application No. 60/214,886, filed on 28-Jun-2000; Application No. 60/217,487, filed on 11-Jul-2000; Application No. 60/225,758, filed on 14-Aug-2000; Application No. 60/220,963, filed on 26-Jul-2000; Application No. 60/217,496, filed on 11-Jul-2000; Application No. 60/225,447, filed on 14-Aug-2000; Application No. 60/218,290, filed on 14-Jul-2000; Application No. 60/225,757, filed on 14-Aug-2000; Application No. 60/226,868, filed on 22-Aug-2000; Application No. 60/216,647, filed on 07-Jul-2000; Application No. 60/225,267, filed on 14-Aug-2000; Application No. 60/216,880, filed on 07-Jul-2000; Application No. 60/225,270, filed on 14-Aug-2000; Application No. 60/251,869, filed on 08-Dec-2000; Application No. 60/235,834, filed on 27-Sep-2000; Application No. 60/234,274, filed on 21-Sep-2000; Application No. 60/234,223, filed on 21-Sep-2000; Application No. 60/228,924, filed on 30-Aug-2000; Application No. 60/224,518, filed on 14-Aug-2000; Application No. 60/236,369, filed on 29-Sep-2000; Application No. 60/224,519, filed on 14-Aug-2000; Application No. 60/220,964, filed on 26-Jul-2000; Application No. 60/241,809, filed on 20-Oct-2000; Application No. 60/249,299, filed on 17-Nov-2000; Application No. 60/236,327, filed on 29-Sep-2000; Application No. 60/241,785, filed on 20-Oct-2000; Application No. 60/244,617, filed on 01-Nov-2000; Application No. 60/225,268, filed on 14-Aug-2000; Application No. 60/236,368, filed on 29-Sep-2000; Application No. 60/251,856, filed on 08-Dec-2000; Application No. 60/251,868, filed on 08-Dec-2000; Application No.

60/229,344, filed on 01-Sep-2000; Application No. 60/234,997, filed on 25-Sep-2000; Application No. 60/229,343, filed on 01-Sep-2000; Application No. 60/229,345, filed on 01-Sep-2000; Application No. 60/229,287, filed on 01-Sep-2000; Application No. 60/229,513, filed on 05-Sep-2000; Application No. 60/231,413, filed on 08-Sep-2000; Application No. 60/229,509, filed on 05-Sep-2000; Application No. 60/236,367, filed on 29-Sep-2000; Application No. 60/237,039, filed on 02-Oct-2000; Application No. 60/237,038, filed on 02-Oct-2000; Application No. 60/236,370, filed on 29-Sep-2000; Application No. 60/236,802, filed on 02-Oct-2000; Application No. 60/237,037, filed on 02-Oct-2000; Application No. 60/237,040, filed on 02-Oct-2000; Application No. 60/240,960, filed on 20-Oct-2000; Application No. 60/239,935, filed on 13-Oct-2000; Application No. 60/239,937, filed on 13-Oct-2000; Application No. 60/241,787, filed on 20-Oct-2000; Application No. 60/246,474, filed on 08-Nov-2000; Application No. 60/246,532, filed on 08-Nov-2000; Application No. 60/249,216, filed on 17-Nov-2000; Application No. 60/249,210, filed on 17-Nov-2000; Application No. 60/226,681, filed on 22-Aug-2000; Application No. 60/225,759, filed on 14-Aug-2000; Application No. 60/225,213, filed on 14-Aug-2000; Application No. 60/227,182, filed on 22-Aug-2000; Application No. 60/225,214, filed on 14-Aug-2000; Application No. 60/235,836, filed on 27-Sep-2000; Application No. 60/230,438, filed on 06-Sep-2000; Application No. 60/215,135, filed on 30-Jun-2000; Application No. 60/225,266, filed on 14-Aug-2000; Application No. 60/249,218, filed on 17-Nov-2000; Application No. 60/249,208, filed on 17-Nov-2000; Application No. 60/249,213, filed on 17-Nov-2000; Application No. 60/249,212, filed on 17-Nov-2000; Application No. 60/249,207, filed on 17-Nov-2000; Application No. 60/249,245, filed on 17-Nov-2000; Application No. 60/249,244, filed on 17-Nov-2000; Application No. 60/249,217, filed on 17-Nov-2000; Application No. 60/249,211, filed on 17-Nov-2000; Application No. 60/249,215, filed on 17-Nov-2000; Application No. 60/249,264, filed on 17-Nov-2000; Application No. 60/249,214, filed on 17-Nov-2000; Application No. 60/249,297, filed on 17-Nov-2000; Application No. 60/232,400, filed on 14-Sep-2000; Application No. 60/231,242, filed on 08-Sep-2000; Application No. 60/232,081, filed on 08-Sep-2000; Application No. 60/232,080, filed on 08-Sep-2000; Application No. 60/231,414, filed on 08-Sep-2000; Application No. 60/231,244, filed on 08-Sep-2000; Application No. 60/233,064, filed on 14-Sep-2000; Application No. 60/233,063, filed on 14-Sep-2000; Application No. 60/232,397, filed on

14-Sep-2000; Application No. 60/232,399, filed on 14-Sep-2000; Application No. 60/232,401, filed on 14-Sep-2000; Application No. 60/241,808, filed on 20-Oct-2000; Application No. 60/241,826, filed on 20-Oct-2000; Application No. 60/241,786, filed on 20-Oct-2000; Application No. 60/241,221, filed on 20-Oct-2000; Application No. 60/246,475, filed on 08-Nov-2000; Application No. 60/231,243, filed on 08-Sep-2000; Application No. 60/233,065, filed on 14-Sep-2000; Application No. 60/232,398, filed on 14-Sep-2000; Application No. 60/234,998, filed on 25-Sep-2000; Application No. 60/246,477, filed on 08-Nov-2000; Application No. 60/246,528, filed on 08-Nov-2000; Application No. 60/246,525, filed on 08-Nov-2000; Application No. 60/246,476, filed on 08-Nov-2000; Application No. 60/246,526, filed on 08-Nov-2000; Application No. PT172, filed on 17-Nov-2000; Application No. 60/246,527, filed on 08-Nov-2000; Application No. 60/246,523, filed on 08-Nov-2000; Application No. 60/246,524, filed on 08-Nov-2000; Application No. 60/246,478, filed on 08-Nov-2000; Application No. 60/246,609, filed on 08-Nov-2000; Application No. 60/246,613, filed on 08-Nov-2000; Application No. 60/249,300, filed on 17-Nov-2000; Application No. 60/249,265, filed on 17-Nov-2000; Application No. 60/246,610, filed on 08-Nov-2000; Application No. 60/246,611, filed on 08-Nov-2000; Application No. 60/230,437, filed on 06-Sep-2000; Application No. 60/251,990, filed on 08-Dec-2000; Application No. 60/251,988, filed on 05-Dec-2000; Application No. 60/251,030, filed on 05-Dec-2000; Application No. 60/251,479, filed on 06-Dec-2000; Application No. PJ005, filed on 05-Dec-2000; Application No. PJ006, filed on 01-Dec-2000; Application No. 60/251,989, filed on 08-Dec-2000; Application No. 60/250,391, filed on 01-Dec-2000; and Application No. 60/254,097, filed on 11-Dec-2000.

**[01281]** Moreover, the microfiche copy and the corresponding computer readable form of the Sequence Listing of U.S. Application Serial No. 60/179,065, and the hard copy of and the corresponding computer readable form of the Sequence Listing of U.S. Application Serial No. 60/180,628 are also incorporated herein by reference in their entireties.

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)

10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209059

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐


**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).  
Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		



**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209060

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209061

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer <div style="text-align: center; font-family: cursive;">ALT</div>			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209062

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)

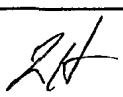
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (including postal code and country)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209063

**C. ADDITIONAL INDICATIONS** (leave blank if not applicable)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (if the indications are not for all designated States)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (leave blank if not applicable)

The indications listed below will be submitted to the international Bureau later (specify the general nature of the indications e.g., "Accession Number of Deposit")

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (including postal code and country)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209064

**C. ADDITIONAL INDICATIONS** (leave blank if not applicable)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (if the indications are not for all designated States)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (leave blank if not applicable)

The indications listed below will be submitted to the international Bureau later (specify the general nature of the indications e.g., "Accession Number of Deposit")

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)

10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209065

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐


**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC). Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209066

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)

10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209067

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)

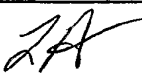
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		



**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209068

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

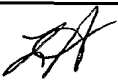
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (including postal code and country)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

May 20, 1997

Accession Number

209069

**C. ADDITIONAL INDICATIONS** (leave blank if not applicable)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (if the indications are not for all designated States)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (leave blank if not applicable)

The indications listed below will be submitted to the international Bureau later (specify the general nature of the indications e.g., "Accession Number of Deposit")

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 12, 1998

Accession Number

209579

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)

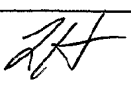
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 12, 1998

Accession Number

209578

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

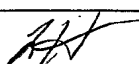
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (*including postal code and country*)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

July 16, 1998

Accession Number

203067

**C. ADDITIONAL INDICATIONS** (*leave blank if not applicable*)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (*if the indications are not for all designated States*)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (*leave blank if not applicable*)

The indications listed below will be submitted to the international Bureau later (*specify the general nature of the indications e.g., "Accession Number of Deposit"*)

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

July 16, 1998

Accession Number

203068

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/>	This sheet was received with the international application			<input type="checkbox"/>	This sheet was received by the International Bureau on:
Authorized officer			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

February 1, 1999

Accession Number

203609

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

February 1, 1999

Accession Number

203610

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

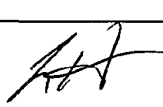
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		



**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*

10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

November 17, 1998

Accession Number

203485

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

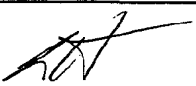
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

June 18, 1999

Accession Number

PTA-252

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

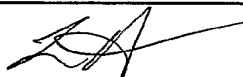
For receiving Office use only

For International Bureau use only

☒ This sheet was received with the international application

☐ This sheet was received by the International Bureau on:

Authorized officer



Authorized officer

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*

10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

June 18, 1999

Accession Number

PTA-253

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

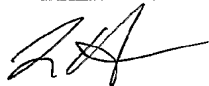
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at Table 6.

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

December 22, 1999

Accession Number

PTA-1081

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐


**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC). Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

October 5, 2000

Accession Number

PTA-2574

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

For receiving Office use only

For International Bureau use only

☒ This sheet was received with the international application

☐ This sheet was received by the International Bureau on:

Authorized officer



Authorized officer

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

October 5, 2000

Accession Number

PTA-2575

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

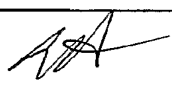
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

**A.** The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 5, 2001

Accession Number

(HGS reference code TS-1)

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

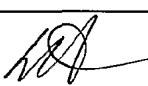
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (including postal code and country)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 5, 2001

Accession Number

(HGS reference code TS-2)

**C. ADDITIONAL INDICATIONS** (leave blank if not applicable)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (if the indications are not for all designated States)


Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (leave blank if not applicable)

The indications listed below will be submitted to the international Bureau later (specify the general nature of the indications e.g., "Accession Number of Deposit")

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		



**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☒

Name of depositary institution: American Type Culture Collection

Address of depositary institution (including postal code and country)  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 5, 2001

Accession Number

(HGS reference code AC-1)

**C. ADDITIONAL INDICATIONS** (leave blank if not applicable)

This information is continued on an additional sheet ☐

**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** (if the indications are not for all designated States)

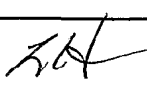
Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC).

Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** (leave blank if not applicable)

The indications listed below will be submitted to the international Bureau later (specify the general nature of the indications e.g., "Accession Number of Deposit")

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**INDICATIONS RELATING TO A DEPOSITED MICROORGANISM  
OR OTHER BIOLOGICAL MATERIAL**

(PCT Rule 13bis)

A. The indications made below relate to the deposited microorganism or other biological material referred to in the description at page 16, paragraph [033].

**B. IDENTIFICATION OF DEPOSIT**

Further deposits are identified on an additional sheet ☐

Name of depositary institution: American Type Culture Collection

Address of depositary institution *(including postal code and country)*  
10801 University Boulevard  
Manassas, Virginia 20110-2209  
United States of America

Date of deposit

January 5, 2001

Accession Number

(HGS reference code AC-2)

**C. ADDITIONAL INDICATIONS** *(leave blank if not applicable)*

This information is continued on an additional sheet ☐

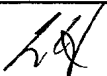
**D. DESIGNATED STATES FOR WHICH INDICATIONS ARE MADE** *(if the indications are not for all designated States)*

Europe

In respect of those designations in which a European Patent is sought a sample of the deposited microorganism will be made available until the publication of the mention of the grant of the European patent or until the date on which the application has been refused or withdrawn or is deemed to be withdrawn, only by the issue of such a sample to an expert nominated by the person requesting the sample (Rule 28(4) EPC). Continued on additional sheets

**E. SEPARATE FURNISHING OF INDICATIONS** *(leave blank if not applicable)*

The indications listed below will be submitted to the international Bureau later *(specify the general nature of the indications e.g., "Accession Number of Deposit")*

	For receiving Office use only			For International Bureau use only	
<input checked="" type="checkbox"/> This sheet was received with the international application			<input type="checkbox"/> This sheet was received by the International Bureau on:		
Authorized officer 			Authorized officer		

**ATCC Deposit No.:** 209059, 209060, 209061, 209062, 209063, 209064, 209065, 209066, 209067, 209068, 209069, 209579, 209578, 203067, 203068, 203609, 203610, 203485, PTA-252, PTA-253, PTA-1081, PTA-2574, PTA-2575, TS-1, TS-2, AC-1, AC-2

#### **CANADA**

The applicant requests that, until either a Canadian patent has been issued on the basis of an application or the application has been refused, or is abandoned and no longer subject to reinstatement, or is withdrawn, the Commissioner of Patents only authorizes the furnishing of a sample of the deposited biological material referred to in the application to an independent expert nominated by the Commissioner, the applicant must, by a written statement, inform the International Bureau accordingly before completion of technical preparations for publication of the international application.

#### **NORWAY**

The applicant hereby requests that the application has been laid open to public inspection (by the Norwegian Patent Office), or has been finally decided upon by the Norwegian Patent Office without having been laid open inspection, the furnishing of a sample shall only be effected to an expert in the art. The request to this effect shall be filed by the applicant with the Norwegian Patent Office not later than at the time when the application is made available to the public under Sections 22 and 33(3) of the Norwegian Patents Act. If such a request has been filed by the applicant, any request made by a third party for the furnishing of a sample shall indicate the expert to be used. That expert may be any person entered on the list of recognized experts drawn up by the Norwegian Patent Office or any person approved by the applicant in the individual case.

#### **AUSTRALIA**

The applicant hereby gives notice that the furnishing of a sample of a microorganism shall only be effected prior to the grant of a patent, or prior to the lapsing, refusal or withdrawal of the application, to a person who is a skilled addressee without an interest in the invention (Regulation 3.25(3) of the Australian Patents Regulations).

#### **FINLAND**

The applicant hereby requests that, until the application has been laid open to public inspection (by the National Board of Patents and Regulations), or has been finally decided upon by the National Board of Patents and Registration without having been laid open to public inspection, the furnishing of a sample shall only be effected to an expert in the art.

**ATCC Deposit No.:** 209059, 209060, 209061, 209062, 209063, 209064, 209065, 209066, 209067, 209068, 209069, 209579, 209578, 203067, 203068, 203609, 203610, 203485, PTA-252, PTA-253, PTA-1081, PTA-2574, PTA-2575, TS-1, TS-2, AC-1, AC-2

## **UNITED KINGDOM**

The applicant hereby requests that the furnishing of a sample of a microorganism shall only be made available to an expert. The request to this effect must be filed by the applicant with the International Bureau before the completion of the technical preparations for the international publication of the application.

## **DENMARK**

The applicant hereby requests that, until the application has been laid open to public inspection (by the Danish Patent Office), or has been finally decided upon by the Danish Patent office without having been laid open to public inspection, the furnishing of a sample shall only be effected to an expert in the art. The request to this effect shall be filed by the applicant with the Danish Patent Office not later than at the time when the application is made available to the public under Sections 22 and 33(3) of the Danish Patents Act. If such a request has been filed by the applicant, any request made by a third party for the furnishing of a sample shall indicate the expert to be used. That expert may be any person entered on a list of recognized experts drawn up by the Danish Patent Office or any person by the applicant in the individual case.

## **SWEDEN**

The applicant hereby requests that, until the application has been laid open to public inspection (by the Swedish Patent Office), or has been finally decided upon by the Swedish Patent Office without having been laid open to public inspection, the furnishing of a sample shall only be effected to an expert in the art. The request to this effect shall be filed by the applicant with the International Bureau before the expiration of 16 months from the priority date (preferably on the Form PCT/RO/134 reproduced in annex Z of Volume I of the PCT Applicant's Guide). If such a request has been filed by the applicant any request made by a third party for the furnishing of a sample shall indicate the expert to be used. That expert may be any person entered on a list of recognized experts drawn up by the Swedish Patent Office or any person approved by a applicant in the individual case.

## **NETHERLANDS**

The applicant hereby requests that until the date of a grant of a Netherlands patent or until the date on which the application is refused or withdrawn or lapsed, the microorganism shall be made available as provided in the 31F(1) of the Patent Rules only by the issue of a sample to an expert. The request to this effect must be furnished by the applicant with the Netherlands Industrial Property Office before the date on which the application is made available to the public under Section 22C or Section 25 of the Patents Act of the Kingdom of the Netherlands, whichever of the two dates occurs earlier.

*What Is Claimed Is:*

1. An isolated nucleic acid molecule comprising a polynucleotide having a nucleotide sequence at least 95% identical to a sequence selected from the group consisting of:
  - (a) a polynucleotide fragment of SEQ ID NO:X or a polynucleotide fragment of the cDNA sequence contained in Clone ID NO:Z, which is hybridizable to SEQ ID NO:X;
  - (b) a polynucleotide encoding a polypeptide fragment of SEQ ID NO:Y or a polypeptide fragment encoded by the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X;
  - (c) a polynucleotide encoding a polypeptide fragment of a polypeptide encoded by SEQ ID NO:X or a polypeptide fragment encoded by the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X;
  - (d) a polynucleotide encoding a polypeptide domain of SEQ ID NO:Y or a polypeptide domain encoded by the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X;
  - (e) a polynucleotide encoding a polypeptide epitope of SEQ ID NO:Y or a polypeptide epitope encoded by the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X;
  - (f) a polynucleotide encoding a polypeptide of SEQ ID NO:Y or the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X, having biological activity;
  - (g) a polynucleotide which is a variant of SEQ ID NO:X;
  - (h) a polynucleotide which is an allelic variant of SEQ ID NO:X;
  - (i) a polynucleotide which encodes a species homologue of the SEQ ID NO:Y;
  - (j) a polynucleotide capable of hybridizing under stringent conditions to any one of the polynucleotides specified in (a)-(i), wherein said polynucleotide does not hybridize under stringent conditions to a nucleic acid molecule having a nucleotide sequence of only A residues or of only T residues.
2. The isolated nucleic acid molecule of claim 1, wherein the polynucleotide fragment comprises a nucleotide sequence encoding a protein.

3. The isolated nucleic acid molecule of claim 1, wherein the polynucleotide fragment comprises a nucleotide sequence encoding the sequence identified as SEQ ID NO:Y or the polypeptide encoded by the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X.

4. The isolated nucleic acid molecule of claim 1, wherein the polynucleotide fragment comprises the entire nucleotide sequence of SEQ ID NO:X or the cDNA sequence contained in cDNA Clone ID NO:Z, which is hybridizable to SEQ ID NO:X.

5. The isolated nucleic acid molecule of claim 2, wherein the nucleotide sequence comprises sequential nucleotide deletions from either the C-terminus or the N-terminus.

6. The isolated nucleic acid molecule of claim 3, wherein the nucleotide sequence comprises sequential nucleotide deletions from either the C-terminus or the N-terminus.

7. A recombinant vector comprising the isolated nucleic acid molecule of claim 1.

8. A method of making a recombinant host cell comprising the isolated nucleic acid molecule of claim 1.

9. A recombinant host cell produced by the method of claim 8.

10. The recombinant host cell of claim 9 comprising vector sequences.

11. An isolated polypeptide comprising an amino acid sequence at least 90% identical to a sequence selected from the group consisting of:

(a) a polypeptide fragment of SEQ ID NO:Y or the encoded sequence contained in cDNA Clone ID NO:Z;

(b) a polypeptide fragment of SEQ ID NO:Y or the encoded sequence contained in cDNA Clone ID NO:Z, having biological activity;

- (c) a polypeptide domain of SEQ ID NO:Y or the encoded sequence contained in cDNA Clone ID NO:Z;
- (d) a polypeptide epitope of SEQ ID NO:Y or the encoded sequence contained in cDNA Clone ID NO:Z;
- (e) a full length protein of SEQ ID NO:Y or the encoded sequence contained in cDNA Clone ID NO:Z;
- (f) a variant of SEQ ID NO:Y;
- (g) an allelic variant of SEQ ID NO:Y; or
- (h) a species homologue of the SEQ ID NO:Y.

12. The isolated polypeptide of claim 11, wherein the full length protein comprises sequential amino acid deletions from either the C-terminus or the N-terminus.
13. An isolated antibody that binds specifically to the isolated polypeptide of claim 11.
14. A recombinant host cell that expresses the isolated polypeptide of claim 11.
15. A method of making an isolated polypeptide comprising:
  - (a) culturing the recombinant host cell of claim 14 under conditions such that said polypeptide is expressed; and
  - (b) recovering said polypeptide.
16. The polypeptide produced by claim 15.
17. A method for preventing, treating, or ameliorating a medical condition, comprising administering to a mammalian subject a therapeutically effective amount of the polynucleotide of claim 1.
18. A method of diagnosing a pathological condition or a susceptibility to a pathological condition in a subject comprising:
  - (a) determining the presence or absence of a mutation in the polynucleotide of claim 1; and

(b) diagnosing a pathological condition or a susceptibility to a pathological condition based on the presence or absence of said mutation.

19. A method of diagnosing a pathological condition or a susceptibility to a pathological condition in a subject comprising:

(a) determining the presence or amount of expression of the polypeptide of claim 11 in a biological sample; and

(b) diagnosing a pathological condition or a susceptibility to a pathological condition based on the presence or amount of expression of the polypeptide.

20. A method for identifying a binding partner to the polypeptide of claim 11 comprising:

(a) contacting the polypeptide of claim 11 with a binding partner; and

(b) determining whether the binding partner effects an activity of the polypeptide.

21. The gene corresponding to the cDNA sequence of SEQ ID NO:Y.

22. A method of identifying an activity in a biological assay, wherein the method comprises:

(a) expressing SEQ ID NO:X in a cell;

(b) isolating the supernatant;

(c) detecting an activity in a biological assay; and

(d) identifying the protein in the supernatant having the activity.

23. The product produced by the method of claim 20.

24. A method for preventing, treating, or ameliorating a medical condition, comprising administering to a mammalian subject a therapeutically effective amount of the polypeptide of claim 11.



(19) World Intellectual Property Organization  
International Bureau(43) International Publication Date  
9 August 2001 (09.08.2001)

PCT

(10) International Publication Number  
**WO 01/57182 A3**

(51) International Patent Classification <sup>7</sup> : C07H 21/02	C12Q 1/68,	60/225,758	14 August 2000 (14.08.2000)	US	
		60/225,267	14 August 2000 (14.08.2000)	US	
		60/225,214	14 August 2000 (14.08.2000)	US	
(21) International Application Number:	PCT/US01/01354	60/226,279	18 August 2000 (18.08.2000)	US	
		60/226,868	22 August 2000 (22.08.2000)	US	
(22) International Filing Date:	17 January 2001 (17.01.2001)	60/227,182	22 August 2000 (22.08.2000)	US	
		60/226,681	22 August 2000 (22.08.2000)	US	
(25) Filing Language:	English	60/227,009	23 August 2000 (23.08.2000)	US	
		60/228,924	30 August 2000 (30.08.2000)	US	
(26) Publication Language:	English	60/229,344	1 September 2000 (01.09.2000)	US	
		60/229,343	1 September 2000 (01.09.2000)	US	
(30) Priority Data:		60/229,287	1 September 2000 (01.09.2000)	US	
		60/229,345	1 September 2000 (01.09.2000)	US	
60/179,065	31 January 2000 (31.01.2000)	US	60/229,513	5 September 2000 (05.09.2000)	US
60/180,628	4 February 2000 (04.02.2000)	US	60/229,509	5 September 2000 (05.09.2000)	US
60/184,664	24 February 2000 (24.02.2000)	US	60/230,438	6 September 2000 (06.09.2000)	US
60/186,350	2 March 2000 (02.03.2000)	US	60/230,437	6 September 2000 (06.09.2000)	US
60/189,874	16 March 2000 (16.03.2000)	US	60/231,413	8 September 2000 (08.09.2000)	US
60/190,076	17 March 2000 (17.03.2000)	US	60/232,080	8 September 2000 (08.09.2000)	US
60/198,123	18 April 2000 (18.04.2000)	US	60/231,414	8 September 2000 (08.09.2000)	US
60/205,515	19 May 2000 (19.05.2000)	US	60/231,244	8 September 2000 (08.09.2000)	US
60/209,467	7 June 2000 (07.06.2000)	US	60/232,081	8 September 2000 (08.09.2000)	US
60/214,886	28 June 2000 (28.06.2000)	US	60/231,242	8 September 2000 (08.09.2000)	US
60/215,135	30 June 2000 (30.06.2000)	US	60/231,243	8 September 2000 (08.09.2000)	US
60/216,647	7 July 2000 (07.07.2000)	US	60/231,968	12 September 2000 (12.09.2000)	US
60/216,880	7 July 2000 (07.07.2000)	US	60/232,401	14 September 2000 (14.09.2000)	US
60/217,487	11 July 2000 (11.07.2000)	US	60/232,399	14 September 2000 (14.09.2000)	US
60/217,496	11 July 2000 (11.07.2000)	US	60/232,400	14 September 2000 (14.09.2000)	US
60/218,290	14 July 2000 (14.07.2000)	US	60/232,397	14 September 2000 (14.09.2000)	US
60/220,963	26 July 2000 (26.07.2000)	US	60/233,063	14 September 2000 (14.09.2000)	US
60/220,964	26 July 2000 (26.07.2000)	US	60/233,064	14 September 2000 (14.09.2000)	US
60/225,757	14 August 2000 (14.08.2000)	US	60/233,065	14 September 2000 (14.09.2000)	US
60/225,270	14 August 2000 (14.08.2000)	US	60/232,398	14 September 2000 (14.09.2000)	US
60/225,447	14 August 2000 (14.08.2000)	US	60/234,223	21 September 2000 (21.09.2000)	US
60/225,266	14 August 2000 (14.08.2000)	US	60/234,274	21 September 2000 (21.09.2000)	US
60/225,213	14 August 2000 (14.08.2000)	US	60/234,997	25 September 2000 (25.09.2000)	US
60/225,759	14 August 2000 (14.08.2000)	US	60/234,998	25 September 2000 (25.09.2000)	US
60/224,519	14 August 2000 (14.08.2000)	US	60/235,484	26 September 2000 (26.09.2000)	US
60/224,518	14 August 2000 (14.08.2000)	US	60/235,834	27 September 2000 (27.09.2000)	US
60/225,268	14 August 2000 (14.08.2000)	US	60/235,836	27 September 2000 (27.09.2000)	US

[Continued on next page]

(54) Title: NUCLEIC ACIDS, PROTEINS, AND ANTIBODIES

(57) **Abstract:** The present invention relates to novel immune/hematopoietic-related polynucleotides and the polypeptides encoded by these polynucleotides herein collectively known as "immune/hematopoietic antigens", and the use of such immune/hematopoietic antigens for detecting immune/hematopoietic-related diseases and/or disorders, particularly the presence of cancer and cancer metastases of cells of hematopoietic origin. More specifically, isolated immune/hematopoietic associated nucleic acid molecules are provided encoding novel immune/hematopoietic associated polypeptides. Novel immune/hematopoietic polypeptides and antibodies that bind to these polypeptides are provided. Also provided are vectors, host cells, and recombinant and synthetic methods for producing human immune/hematopoietic associated polynucleotides and/or polypeptides. The invention further relates to diagnostic and therapeutic methods useful for diagnosing, treating, preventing and/or prognosing disorders related to the immune system or cells and tissues associated with hematopoiesis, including cancers of cells of hematopoietic origin, and therapeutic methods for treating such disorders. The invention further relates to screening methods for identifying agonists and antagonists of polynucleotides and polypeptides of the invention. The present invention further relates to methods and/or compositions for inhibiting the production and function of the polypeptides of the present invention.



60/236,369	29 September 2000 (29.09.2000)	US	60/256,719	5 December 2000 (05.12.2000)	US
60/236,327	29 September 2000 (29.09.2000)	US	60/251,030	5 December 2000 (05.12.2000)	US
60/236,370	29 September 2000 (29.09.2000)	US	60/251,988	5 December 2000 (05.12.2000)	US
60/236,368	29 September 2000 (29.09.2000)	US	60/251,479	6 December 2000 (06.12.2000)	US
60/236,367	29 September 2000 (29.09.2000)	US	60/251,869	8 December 2000 (08.12.2000)	US
60/237,039	2 October 2000 (02.10.2000)	US	60/251,856	8 December 2000 (08.12.2000)	US
60/237,038	2 October 2000 (02.10.2000)	US	60/251,868	8 December 2000 (08.12.2000)	US
60/237,040	2 October 2000 (02.10.2000)	US	60/251,990	8 December 2000 (08.12.2000)	US
60/237,037	2 October 2000 (02.10.2000)	US	60/251,989	8 December 2000 (08.12.2000)	US
60/236,802	2 October 2000 (02.10.2000)	US	60/254,097	11 December 2000 (11.12.2000)	US
60/239,937	13 October 2000 (13.10.2000)	US	60/259,678	5 January 2001 (05.01.2001)	US
60/239,935	13 October 2000 (13.10.2000)	US			
60/241,785	20 October 2000 (20.10.2000)	US	<b>(71) Applicant (for all designated States except US): HUMAN GENOME SCIENCES, INC.</b> [US/US]; 9410 Key West Avenue, Rockville, MD 20850 (US).		
60/241,809	20 October 2000 (20.10.2000)	US	<b>(72) Inventors; and</b>		
60/240,960	20 October 2000 (20.10.2000)	US	<b>(75) Inventors/Applicants (for US only): ROSEN, Craig, A.</b> [US/US]; 22400 Rolling Hill Lane, Laytonsville, MD 20882 (US). <b>BARASH, Steven, C.</b> [US/US]; 111 Watkins Pond Boulevard #301, Rockville, MD 20850 (US). <b>RUBEN, Steven, M.</b> [US/US]; 18528 Heritage Hills Drive, Olney, MD 20832 (US).		
60/241,787	20 October 2000 (20.10.2000)	US	<b>(74) Agents: HOOVER, Kenley, K.</b> et al.; Human Genome Sciences, Inc., 9410 Key West Avenue, Rockville, MD 20850 (US).		
60/241,808	20 October 2000 (20.10.2000)	US	<b>(81) Designated States (national):</b> AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CR, CU, CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW.		
60/241,221	20 October 2000 (20.10.2000)	US	<b>(84) Designated States (regional):</b> ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).		
60/241,786	20 October 2000 (20.10.2000)	US	<b>Published:</b>		
60/241,826	20 October 2000 (20.10.2000)	US	— with international search report		
60/244,617	1 November 2000 (01.11.2000)	US	— with sequence listing part of description published separately in electronic form and available upon request from the International Bureau		
60/246,474	8 November 2000 (08.11.2000)	US	<b>(88) Date of publication of the international search report:</b>		
60/246,532	8 November 2000 (08.11.2000)	US	28 March 2002		
60/246,476	8 November 2000 (08.11.2000)	US	<i>For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.</i>		
60/246,526	8 November 2000 (08.11.2000)	US			
60/246,475	8 November 2000 (08.11.2000)	US			
60/246,525	8 November 2000 (08.11.2000)	US			
60/246,528	8 November 2000 (08.11.2000)	US			
60/246,527	8 November 2000 (08.11.2000)	US			
60/246,477	8 November 2000 (08.11.2000)	US			
60/246,611	8 November 2000 (08.11.2000)	US			
60/246,610	8 November 2000 (08.11.2000)	US			
60/246,613	8 November 2000 (08.11.2000)	US			
60/246,609	8 November 2000 (08.11.2000)	US			
60/246,478	8 November 2000 (08.11.2000)	US			
60/246,524	8 November 2000 (08.11.2000)	US			
60/246,523	8 November 2000 (08.11.2000)	US			
60/249,299	17 November 2000 (17.11.2000)	US			
60/249,210	17 November 2000 (17.11.2000)	US			
60/249,216	17 November 2000 (17.11.2000)	US			
60/249,217	17 November 2000 (17.11.2000)	US			
60/249,211	17 November 2000 (17.11.2000)	US			
60/249,215	17 November 2000 (17.11.2000)	US			
60/249,218	17 November 2000 (17.11.2000)	US			
60/249,208	17 November 2000 (17.11.2000)	US			
60/249,213	17 November 2000 (17.11.2000)	US			
60/249,212	17 November 2000 (17.11.2000)	US			
60/249,207	17 November 2000 (17.11.2000)	US			
60/249,245	17 November 2000 (17.11.2000)	US			
60/249,244	17 November 2000 (17.11.2000)	US			
60/249,297	17 November 2000 (17.11.2000)	US			
60/249,214	17 November 2000 (17.11.2000)	US			
60/249,264	17 November 2000 (17.11.2000)	US			
60/249,209	17 November 2000 (17.11.2000)	US			
60/249,300	17 November 2000 (17.11.2000)	US			
60/249,265	17 November 2000 (17.11.2000)	US			
60/250,391	1 December 2000 (01.12.2000)	US			
60/250,160	1 December 2000 (01.12.2000)	US			

## INTERNATIONAL SEARCH REPORT

International application No.

PCT/US01/01354

**A. CLASSIFICATION OF SUBJECT MATTER**

IPC(7) : C12Q 1/68; C07H 21/02

US CL : 435/6; 536/23.1

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 435/6; 536/23.1

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)  
STIC Biotech Sequence Homology Search (Commercial and issued-Patent databases)**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	Database EST STIC Biotech Search. Accession Number AW075920.1 NCI-CGAP, October 13, 1999.	1-4, 11-12, 14-16
X	Database PIR_66 STIC Biotech Search. Accession Number A496780 Barford et al. 1993.	1-4, 11-12, 14-16
Y	WATSON et al. The Science Used in the Recombination DNA Industry. Recombinant DNA; A Short Course. New York. W.H. Freeman and Company. 1983, pages 231-241, Chapter 18, see entire document.	5-10, 14-16

☐ Further documents are listed in the continuation of Box C.☐ See patent family annex.

\* Special categories of cited documents:

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier application or patent published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T"

later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X"

document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y"

document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art

"&amp;"

document member of the same patent family

Date of the actual completion of the international search

23 April 2001 (23.04.2001)

Date of mailing of the international search report

03 AUG 2001

Name and mailing address of the ISA/US

Commissioner of Patents and Trademarks

Box PCT

Washington, D.C. 20231

Facsimile No. (703)305-3230

Authorized officer

Young J. Kim

Telephone No. (703) 308-0196

## INTERNATIONAL SEARCH REPORT

International application No.

PCT/US01/01354

**BOX II. OBSERVATIONS WHERE UNITY OF INVENTION IS LACKING** This application contains the following inventions or groups of inventions which are not so linked as to form a single general inventive concept under PCT Rule 13.1. In order for all inventions to be examined, the appropriate additional examination fees must be paid.

Group I, claim(s) 1-12, 14-16, and 21, drawn to cDNA, polypeptides, genes, a method of using the cDNA to make host cells comprising the cDNA, and a method of making the polypeptide.

Group II, claim(s) 13, drawn to an antibody specific for the polypeptides of Group I.

Group III, claim(s) 17, drawn to a therapeutic method of using the cDNA or the polypeptide of Group I.

Group IV, claim(s) 18 and 19, drawn to a diagnostic method of using the cDNA or polypeptide of Group I.

Group V, claim(s) 20, drawn to a method of using the polypeptide of Group I to isolate a binding partner.

Group VI, claim(s) 22, drawn to a method of using the cDNA of Group I to identify the activity of the polypeptide encoded by the cDNA.

Group VII, claim(s) 23, drawn to the binding partner made by the method of Group V.

The inventions listed as Groups I-VII do not relate to a single general inventive concept under PCT Rule 13.1 because, under PCT Rule 13.2, they lack the same or corresponding special technical features for the following reasons: PCT Rule 13.1 and Annex B do not provide for unity of invention between two or more different products or methods of use that share a special technical feature.

In addition, each Group detailed above reads on distinct Group drawn to multiple SEQ ID Numbers. The sequences are distinct because they are unrelated sequences, and a further lack of unity is applied to each Group. The lack of unity is partially waived and the Applicants must further elect 10 SEQ ID Numbers for examination in the elected Group detailed above.

# INTERNATIONAL SEARCH REPORT

International application No.

PCT/US01/01354

## Box I Observations where certain claims were found unsearchable (Continuation of Item 1 of first sheet)

This international report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:

1. ☐ Claim Nos.:  
because they relate to subject matter not required to be searched by this Authority, namely:
2. ☐ Claim Nos.:  
because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out, specifically:
3. ☐ Claim Nos.:  
because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).

## Box II Observations where unity of invention is lacking (Continuation of Item 2 of first sheet)

This International Searching Authority found multiple inventions in this international application, as follows:  
Please See Continuation Sheet

1. ☐ As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims.
2. ☐ As all searchable claims could be searched without effort justifying an additional fee, this Authority did not invite payment of any additional fee.
3. ☐ As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims for which fees were paid, specifically claims Nos.:
4. ☒ No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.: 1-12, 14-16, 21, and SEQ ID Numbers 11 and 9763

Remark on Protest ☐ The additional search fees were accompanied by the applicant's protest.  
☐ No protest accompanied the payment of additional search fees.